

KERRY

Biotechnologies technical manual

**Technologies for cell culture,
fermentation, and diagnostic media**



Introduction

For more than 80 years, we have earned our reputation for reliability and excellence in serving the biotech, pharmaceutical and fermentation markets by delivering unsurpassed product value, enhanced opportunities for supply chain consolidation, and new technology innovation.

We bring together our superior products with new innovative solutions and market-driving alternatives to help our customers succeed in today's challenging global marketplace.

Customers benefit from our application expertise and vast global resources of a \$6 billion global leader. We have the worldwide resources and global technical platform to deliver consistent, high quality products backed by unparalleled service, technical support and formulation customization capabilities.

At Kerry, we make it our business to exceed your expectations every day. We have strong recognition in the Pharma and Biotech sectors, with a history of scientific references totaling more than 500. Our core technical expertise in bio-technology is combined with continued investment in validation, growth performance, regulatory compliance and application support; all of which enable us to meet the requirements of the pharmaceutical industry. Whether your goal is disease treatment, disease prevention, novel fermentation derivatives, or media innovation, Kerry understands the requirements and is the partner you can depend on. We prove it everyday. Let us prove it to you.



CONTENTS

- 4** Manufacturing processes for protein hydrolysates
 > Protein hydrolysates
- 7** Manufacturing processes for yeast specialties
- 8** Applications
 > Cell culture media
 > Vaccines
 > Diagnostic media
 > Pharma fermentation media
 > Classical fermentation media
 > Product evolution
- 14** Bio-science product selection
 > Non-animal typical values
 > Dairy & animal typical values
- 30** Appendix
- 31** List of microorganisms
- 38** Cell line references
- 41** References
- 47** Literature references cited
-

Manufacturing processes for protein hydrolysates

Manufacturing centre

Kerry protein hydrolysates are manufactured primarily in Norwich (NY), but also Rochester (MN) USA and Utrecht, The Netherlands. The existence of these three sites offers unique flexibility in capabilities and contingency. Our sites are supported by product development groups and state of the art pilot plants to accommodate new development, process optimization and customer-specific projects.

Production standards

Kerry is committed to exceeding the cGMP expectations of our pharmaceutical customer base. We set the highest standards and expectations for our systems, written standard operating procedures (SOPs), and documentation.

Our Sheffield™ Brand protein hydrolysates have an extensive quality control program. Testing is not only done on the final product, but also during the manufacturing process. An overview of our standard testing methods can be found on page 30. In addition to the analysis of chemical and physical parameters, the performance of the hydrolysates in application is included in our quality control program.

Kerry has a dedicated Regulatory Group to ensure product compliance to global regulatory requirements and effective reaction to global developments. Regulatory capabilities include:

Regulatory compliance and documentation

- > APHIS inspected and certified
- > Sourcing from countries known to be BSE free
- > No specified risk materials (SRM) used in a manufacturing facility free of specified risk material (SRM)
- > Compliance of milk-derived products to CPMP requirements
- > Formal cleaning validation
- > Availability of certificate of suitability for meat-derived peptones
- > Availability of certificate of origin for all products
- > Traceability of raw materials

Supplier documentation

- > Raw materials from approved suppliers
- > Supplier program
- > Kerry required health certificates for animal-derived materials

Customer care support can also assist with:

- > Confidentiality agreements
- > Import/export documents
- > Completion of questionnaires
- > Customer audits of the manufacturing facilities



Protein hydrolysates

Protein hydrolysates are obtained by (enzymatic or acid) hydrolysis of protein. In this process, proteins are digested into smaller fragments via cleavage of the peptide bond (Fig. 1). Protein hydrolysates are also often referred to as peptides or peptones.

Kerry has many years of experience in protein hydrolysis and has built up in-depth knowledge on the control of the hydrolysis process. When combined with our international expertise in fermentation technology and microbiology, this knowledge gives Kerry the ability to produce protein hydrolysates that are optimally targeted for specific applications, providing unique solutions to our valued customers.

The manufacturing of protein hydrolysates involves several strictly controlled processes, including careful selection of the protein raw materials and enzymes (Fig 2).

The properties of a protein hydrolysate are largely determined by:

1. Protein raw material(s) selection
2. Type and amount of enzyme(s)
3. Reaction conditions during hydrolysis
4. Downstream processing after hydrolysis

Typical protein raw materials for the protein hydrolysates include:

- > Milk
- > Casein
- > Whey protein
- > Animal derived proteins
- > Meat
- > Collagen
- > Vegetable proteins
- > Wheat gluten
- > Soy protein
- > Rice protein
- > Pea protein
- > Cotton seed protein

FIG. 1: THE PEPTIDE BOND


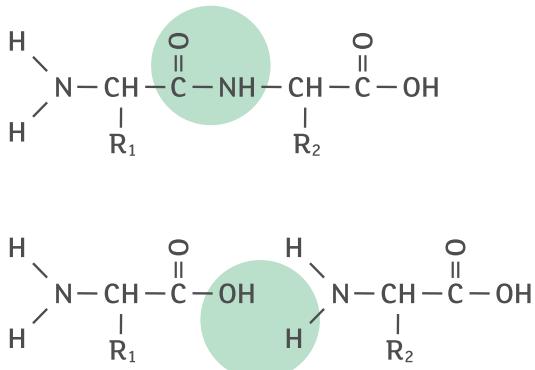


FIG. 2: THE PROTEIN HYDROLYSIS PROCESS

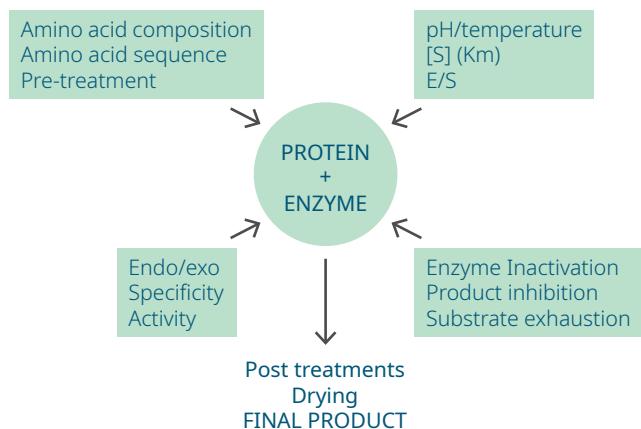

Manufacturing processes for protein hydrolysates *continued*

Protein hydrolysates

For the hydrolysis of these proteins, both endo- and exopeptidases can be used. Endopeptidases cleave internal bonds within the protein chain, and therefore generate (poly) peptides. Exopeptidases hydrolyze proteins and polypeptide chains from the ends of the chain and therefore release predominately amino acids.

The degree of hydrolysis (DH) is a key parameter in the hydrolysis process. It is defined as the percentage of the peptide bonds cleaved, and can be estimated from the ratio between the amino nitrogen (AN) and the total nitrogen (TN). The concentration of free amino acids is also an important parameter. The DH is affected by temperature, pH, time and the concentrations of protein and enzyme. Fig. 3 shows the relationships between the process time, the enzyme-protein ratio and the degree of hydrolysis. Longer hydrolysis and the use of more enzyme result in a higher DH. In-

house enzyme research give Kerry the detailed knowledge of these relationships, and allows for tight control of the hydrolysis process. This enables us to make consistent, high quality products to meet our customers' requirements.

After the hydrolysis step, the enzyme is inactivated, and the hydrolysate may be further processed and purified. Industrially applied processes include centrifugation, plate and frame filtration and membrane filtration (Fig. 4). Centrifugation and plate and frame filtration are used to remove insoluble components to obtain a final product with very high clarity (important for diagnostic media) and filterability (important when the product has to be sterile filtered). Membrane filtration can be utilized to remove larger peptides and endotoxins (important for cell tissue culture media). The final process typically consists of drying the hydrolysate. Kerry, however, can also supply concentrated forms of some hydrolysates.

FIG. 3: THE DEGREE OF HYDROLYSIS IS A KEY PARAMETER IN THE HYDROLYSIS PROCESS


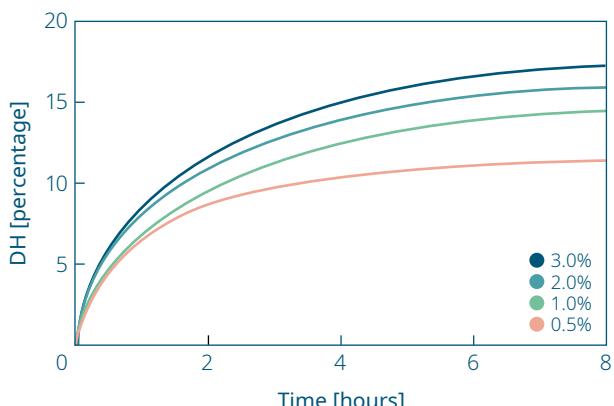


FIG. 4: DOWNSTREAM PROCESSING OF PROTEIN HYDROLYSATES

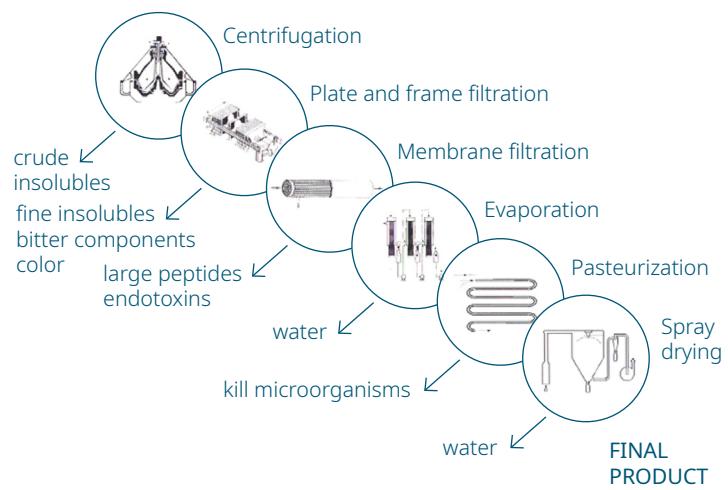

Manufacturing processes of yeast for biotech applications

Manufacturing centre

Kerry's yeast specialties and yeast autolysates are derived from primary grown bakers yeast, propagated aerobically under controlled, reproducible conditions in our yeast fermentation plant in Menstrie, Scotland (UK) and Norwich, NY (USA).

Production standards

For more than thirty years, our factory in Menstrie has successfully supplied pharma customers. The standards demanded by our customer base are exceeded through investment in state-of-the-art hygienically designed equipment and control systems. Cleaning regimes have been developed as part of product safety risk assessments (HACCP) and GMP practices have been put in place throughout the manufacturing processes.

The management of GMP and HACCP is driven through strict auditing, monitoring and validation of CCP's and, where demanded, documented cleaning procedures. The quality system has been in place site-wide since the early 1990's and has provided a structure which delivers continuous improvement to all aspects of the Menstrie operation. Currently, the facility is certified to ISO 9001 which encompasses all activities relating to development, manufacturing, supply and delivery.

Our manufacturing facility in Menstrie uses only animal-free raw materials and provides a Kosher Pareve and Halal production range. This certification is achieved through compliance with the leading global assessment bodies.

Production details

Kerry yeast autolysates and extracts ensure optimal performance in the growth of a broad range of industrially significant microorganisms by providing essential peptides, amino acids, vitamins and minerals, as well as a complex mixture of yeast-derived metabolites. The extracts and autolysates from primary grown yeast (sold by Kerry) offer a distinct performance and consistency advantage over secondary, or brewer's yeast extracts. Brewer's yeast extracts (non Kerry) are produced from the waste streams of the brewing industry and thus are derived from a wide variety of yeast strains and fermentation processes.

After propagation, the baker's yeast cells are harvested and washed via centrifugation. The yeast obtained is then further processed into yeast specialties. In the production of autolyzed yeast, the yeast is heated to induce autolysis, pasteurized and roller dried. The dried product can be packaged directly after sieving as a granular/flake type material or milled to obtain a powder. In the production of yeast extracts, autolysis or digestion is followed by the removal of the insoluble components to give a soluble extract; this is concentrated by vacuum evaporation and filtered (where applicable). Products can be supplied as liquid, paste or non-dusting powder, depending on customer requirements.

Applications

People have used micro-organisms for centuries, e.g. to make beer, wine, cheese and bread. In the last century, many new applications have been added, resulting in a broad-based, ever-growing biotechnology industry. Kerry has over 65 years of experience in peptone optimization for a large range of bacteria, yeasts and molds. We have the knowledge and experience to meet the nutrient requirements of the biotechnology industry, and our products have a reputation for consistency and performance in:

- > cell culture media
- > vaccine media
- > pharma fermentation media
- > diagnostic media
- > classical fermentation media
- > stem cell, regenerative medicine

Cell culture media

Protein hydrolysates are routinely employed to enhance the overall performance of many biopharmaceutical production systems, providing a number of different benefits. Cell viability, cell proliferation and target protein production all may be improved. However, these effects might not be observed concurrently in a given system. As every cell culture platform is unique, so are the protein hydrolysates used to supplement the culture medium. Each may have its own distinctive effect on the growth and productivity of transfected cells. The contribution of protein hydrolysates to the system performance can be influenced by a number of factors including the specific cell line being employed, the raw material used to manufacture the hydrolysate, the manufacturing process, the hydrolysate dosage, and the composition of the basal growth medium.

Supplementation of cell culture media with protein hydrolysates can yield a number of benefits. However, due to the complexity of culture dynamics in supplemented systems, a certain degree of thoughtful experimentation is required to determine the appropriate supplementation scheme that will generate the desired balance among growth, cell viability and target protein titer. Once achieved, the proper supplementation regime can significantly enhance the productivity of any given system.

To develop a better understanding of our cell culture products and how they function in a variety of cell culture systems, we opened the Cell Nutrition Innovation Laboratory. With state-of-the-art equipment and in-house expertise, the laboratory is evidence of our growing capabilities, and the increasing scope of our activities. More importantly, it gives us increased capacity to collaborate with customers on medium optimization and supplement customization to provide process and cell line specific solutions.

© Dennis Kunkel, Microscopy, Inc.

Vaccines

Since 1798 when Dr. David Edward Jenner developed the first vaccine against smallpox, the vaccine industry has gone through several ups and downs. Currently, the vaccine market is enjoying an explosive growth due to development of new technologies, increased awareness, expanded markets, and less government restrictions.

Vaccines can be categorized into human and veterinary. They are classified by how they are produced. The FDA classifies vaccines (prophylactic or therapeutic) based on the sources they come from. These sources include animal, viral, or cellular (microbial or animal cells). Animal sources include fertilized chicken eggs. Virus sourced vaccines are grown in cells or tissue culture.


Kerry offers products that can be used in various stages of vaccine production.

These include:

1. Antigen generation by growing bacteria in fermentation, growing viruses on primary or continuous cell lines, or expression of recombinant proteins derived from disease causing bacteria or viruses in yeast, bacteria, or cell culture. The most common primary cell line used is Chicken Embryonic Fibroblast (CEF) cell line. Human diploid cell lines such as MRC-5 and WI-38 are also widely used. Transformed cell lines including VERO, MDCK, and BHK are also employed. The most frequently used microbes for recombinant bacterial vaccines are E. coli and S. cerevisiae.
2. Antigen isolation and purification
3. Vaccine formulation, which includes addition of adjuvant, stabilizers, and other excipients.

Kerry provides:

- > Protein hydrolysates and yeast extracts for growing microbial, mammalian, and insect cells
- > Animal-component-free (ACF) supplements for reduction or elimination of serum
- > Protein hydrolysates as stabilizers and diluents in vaccine formulation
- > Emulsifiers and surfactants for veterinary vaccine formulation
- > Pharma grade multi-compendial lactose for vaccine stabilization

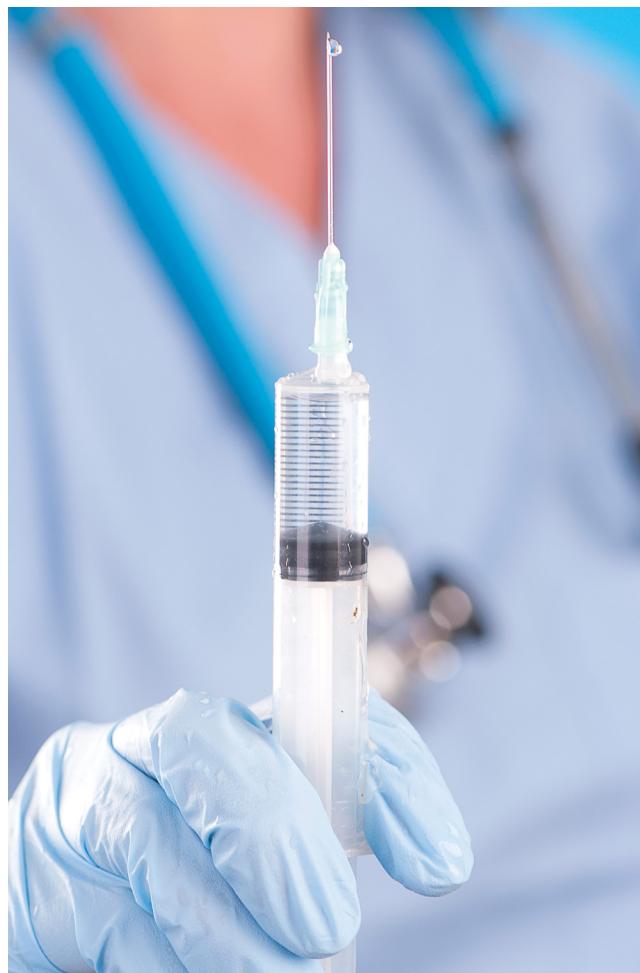

Applications *continued*

Diagnostic medium

Kerry started producing Sheffield™ brand casein hydrolysates for use in laboratory media in 1940. Since then, many of our high quality peptones have been used in clinical diagnostics and quality control media formulations.

In 2004 we expanded our capabilities with the introduction of a fully equipped diagnostic application laboratory, located in Beloit, WI. The application lab enables us to specifically select and optimize peptones, yeast extracts and combinations for specific media applications (general purpose, chromogenic and blood agar medium) and allows us to optimize new products before we start commercializing them.

Kerry diagnostic media components: Typical characteristics

- > Sourcing from countries known to be BSE free
- > Low lot-to-lot variability
- > High degree of clarity for accurate colony identification
- > Unique content of peptides, growth factors, trace elements
- > Minimal color contribution to medium
- > Free of cytotoxic or hemolytic components
- > High solubility
- > Near neutral pH values
- > Promote steady ionic environment
- > Heat stable during autoclaving
- > Shelf-life stability



Pharma fermentation medium

Kerry has a long history of being the number 1 global supplier of nitrogen sources into medium for microbial expression of recombinant proteins. In 1977, N-Z-Amine™ A (enzymatically hydrolysed casein) was used in pioneering work with Charon phages as vectors for cloning DNA in *Escherichia coli* (Blattner, 1977). Since that time, countless researchers have reported the use of N-Z-Amine™ A in both bench level and large scale production of pharmaceuticals.

Today we are providing peptones in the production and/or feed medium of 50% of the top 12 blockbuster drugs produced by *Escherichia coli* (data 2009).

Although most commercially available biopharmaceuticals are still manufactured using casein-based peptones, there is a clear trend towards non-animal medium in this market. In the 1940's, Kerry introduced the first non-animal peptones based on soy based raw materials and expanded this product line in the early 90's with yeast extracts and in 2002 with hydrolysates based on pea, rice, wheat gluten and cottonseed raw materials.

In 2005, we opened a GMO application lab where we work with our in-house recombinant model systems and customer strains to optimize their media to reach maximum productivity.

In our application laboratory, we have the resources to initiate and assist customers with media optimization projects. Medium screening studies can be carried out at three levels:

- > Quick screening for growth optimization in microtiter plates
- > Screening for growth and product formation in O₂ and pH-regulated shake flasks
- > Translation of shake flask results to fermentors to ensure customer compatibility


Kerry is the only medium ingredients supplier worldwide that manufactures both yeast extracts and peptones. This provides us with the unique opportunity to supply nitrogen source solutions optimized for maximum expression of recombinant product in specific expression systems.

It is our pleasure to complement your own expertise with ours, in designing a superior culture medium for your microbial expression system.

During registration of your new products or biosimilars, our dedicated global Kerry regulatory team is pleased to assist you with the information needed for your registration files.



Applications *continued*

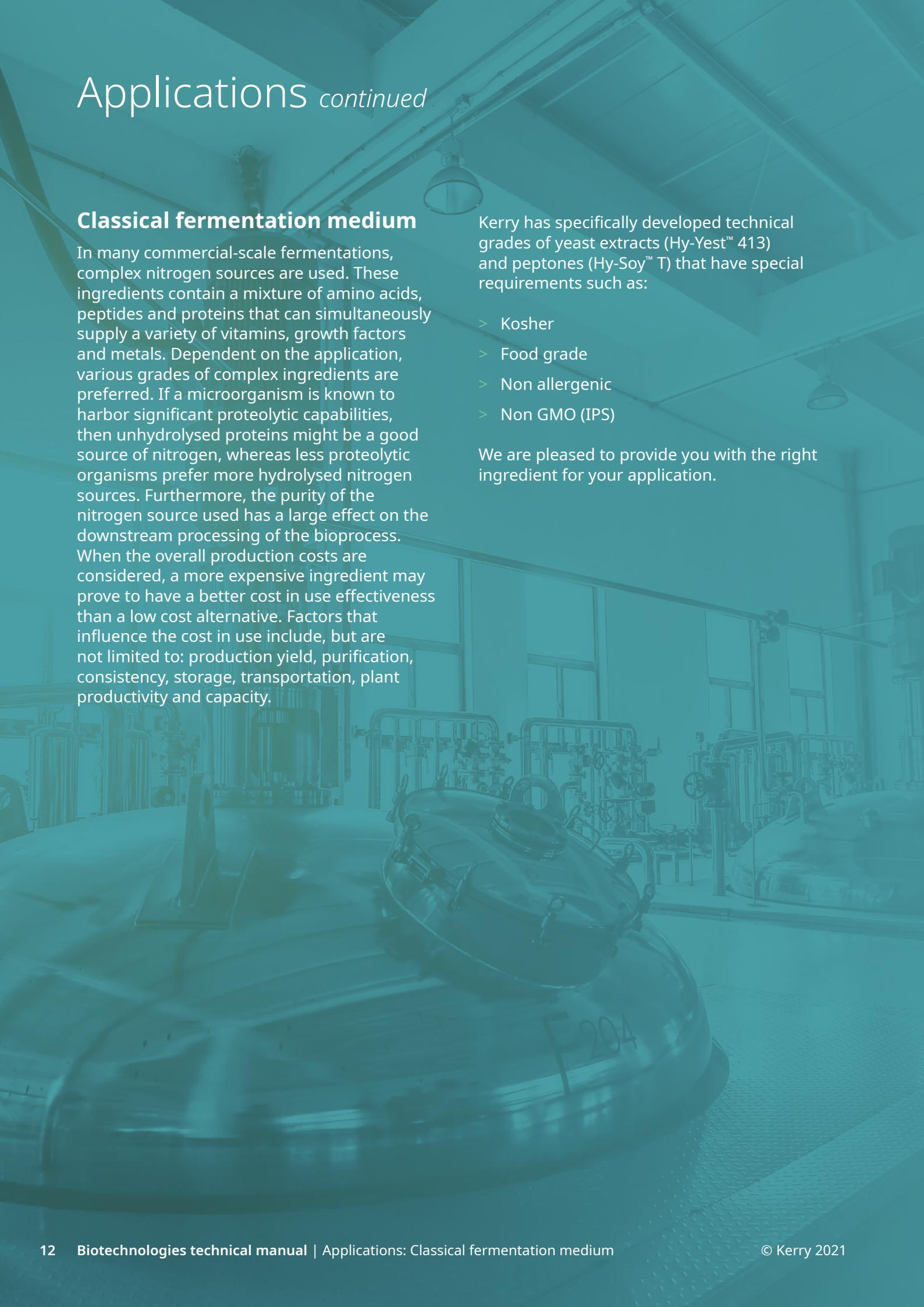A large, stainless steel industrial fermentation tank is shown in the foreground, with the identifier "F204" visible on its side. In the background, there are rows of similar tanks and complex piping systems, all set within a large, modern industrial facility.

Classical fermentation medium

In many commercial-scale fermentations, complex nitrogen sources are used. These ingredients contain a mixture of amino acids, peptides and proteins that can simultaneously supply a variety of vitamins, growth factors and metals. Dependent on the application, various grades of complex ingredients are preferred. If a microorganism is known to harbor significant proteolytic capabilities, then unhydrolysed proteins might be a good source of nitrogen, whereas less proteolytic organisms prefer more hydrolysed nitrogen sources. Furthermore, the purity of the nitrogen source used has a large effect on the downstream processing of the bioprocess. When the overall production costs are considered, a more expensive ingredient may prove to have a better cost in use effectiveness than a low cost alternative. Factors that influence the cost in use include, but are not limited to: production yield, purification, consistency, storage, transportation, plant productivity and capacity.

Kerry has specifically developed technical grades of yeast extracts (Hy-Yest™ 413) and peptones (Hy-Soy™ T) that have special requirements such as:

- > Kosher
- > Food grade
- > Non allergenic
- > Non GMO (IPS)

We are pleased to provide you with the right ingredient for your application.

Product evolution

New generation of hydrolysates

Traditionally, protein hydrolysates have been produced via the proteolytic digestion of a given raw material, followed by iterative filtration and evaporation steps ultimately yielding a spray-dried powder. Certain factors inherent to this type of processing may have limiting or negative effects on the performance of the final product in eukaryotic cell systems. Historical data have shown a direct link between lot-to-lot variability in product composition, and lot-to-lot variability in hydrolysate performance in cell culture systems.

UltraPep™ protein hydrolysates are manufactured using a novel approach to the digestion and processing of plant-derived raw materials that yields high-quality hydrolysates, with increased lot-to-lot consistency and enhanced performance characteristics, for use in eukaryotic biopharmaceutical production systems.

Animal-Component-Free (ACF) defined supplements

Much success has been achieved using plant-derived protein hydrolysates in combination with other supplements such as albumin, normally derived from FBS or other animal sources. To further meet the challenging needs of the market by enhance performance and consistency, Kerry offers high-quality and cost effective recombinant supplements for the supplementation of cell culture media.

Complete supplement solutions

Leveraging our technical and manufacturing expertise of cell culture medium supplements, our application labs developed carefully optimized formulations to further enhance the performance of a variety of cell lines. These highly optimized systems are rapid cost effective supplements that streamline the screening process and while providing customers with requirement options such as protein free or chemically defined.



Bio-science product selection

PRODUCT SELECTION: CELL CULTURE

Supplement System	Cell Lines			
	CHO	CHO Fed Batch	HEK	PerC6
Bio Production Systems				
Sheff-CHO ACF	★	●	★	●
Sheff-CHO Plus ACF	★	●		●
Sheff-CHO PF ACF	★	●		●
Sheff-CHO Plus PF ACF	★	●	★	●
Sheff-Pulse I		★		
Sheff-Pulse II		★		

* Insect cells can be used for bioproduction. We recommend the same products that are used for vaccines in the table below

Supplement System	Cell Lines									
	HEK	VERO	PerC6	MRC-5	BHK	CEF	MDCK	SF-9	High-Five	S2
Vaccines										
Sheff-Vax ACF	●	●	★					★	★	★
Sheff-Vax PF ACF	●	●	★		★	★		★	★	★
Sheff-Vax Plus ACF	●	●	★	●			★	★	★	★
Sheff-Vax Plus PF ACF	●	●	★			★		★	★	★
Sheff-Vax Plus PF ACF VP	★	★	★	★	★	★	★	★	★	★
Sheff-Vax PB10					★	★				

Chart Key: ● Effective ★ Most effective

Notes: Plus: includes supplementation with trace elements

PF: protein free

Product	Cell Lines													
	CHO	CHO Fed Batch	HEK	Hybridoma	VERO	PerC6	MRC-5	BHK	CEF	MDCK	SF-9	High-Five	S2	Stem Cells
Defined Supplements														
Sheffield™ rTrypsin ACF	● *			● *			●	●	●	● *	●	● *	● *	●
Sheffield™ rAlbumin ACF**	★	●	●		●	●	●	●	●	●	●	●	●	●
Sheffield™ rAlbumin EG														●
Complex Inhibitor														
Sheffield™ TrypStop ACF					●	●	●		●	●				●
Ultra-Filtered, Non-Animal Hydrolysates														
HyPep™ 1510	●	●	●		●	●	●	●	●	●	●			
HyPep™ 1511**		●												
HyPep™ 1512	●	●	●		●	●	●	●	●	●	●	●		
HyPep™ 4601N	●		●		●	●			●	●				
HyPep™ 5603N**	●				●	●	●	●			●			
HyPep™ 7504	●	●	●		●	●	●	●	●	●	●			●
UltraPep™ Soy**	●	●	●			●	●			●				
Ultra-Filtered Yeast Extracts														
HyPep™ YE	●			●	●	●	●				●	●	●	●
UltraPep™ YE	●	●		●	●		●				●	●	●	●
Primatone™ HS/UF														
Primatone™ RL/UF								●	●	●	●			

Chart Key: ● Effective ★ Most effective

Notes: CHO: K-1, S, DG44, dhfr-, DXB11

Hybridoma: Myeloma, NSO, SP2/0, NS1

* adherent

** also recommended for Vaccines

Bio-science product selection *continued*

PRODUCT SELECTION: FERMENTATION

We provide total nitrogen sources optimized for your application

Microorganism	Application	Non-Animal										Dairy & Meat																											
		Hy-Express™ System I	Hy-Express™ System II	Hy-Express™ System IV	Hy-Express System V	Hy-Express System VII	Amiso® **	Amiso® BF **	Hy-Express™ Peptone I	Hy-Express™ Peptone II	Hy-Pea® 7404 **	Hy-Pea® 7477 **	Hy-Ric® 5377 **	Hy-Soy™ **	Hy-Soy™ T **	Hy-Soy™ Bl 4	N-Z Soy® K	Sheff TSA/TSB Complete	Hy-Yest® 101 **	Hy-Yest® 412 **	Hy-Yest® 413 **	Hy-Yest® 444 **	Hy-Yest® 455 **	Hy-Yest® 466 **	Hy-Yest® 502 **	Hy-Yest® 503 **	Hy-Yest® 504 **	Hy-Yest® 555 **	Sheftone™ R	Amicase™	Hy-Case™ Amino **	Hy-Case™ Sf	N-Z-Amine™ A	N-Z-Amine™ AS	N-Z-Amine™ EKC **	N-Z-Case TT	Edamin™ F	Pepticas®	Peptonized Milk
<i>Aspergillus</i>	Various																																						
<i>Bacillus</i>	Various	A	⊕ ⊕ ⊕					●		●			●	●	●	●	●			●	●	●	●	●	●	●													
<i>Bifidobacterium sp.</i>	Cultures									●	●									●	●	●																	
<i>Bordetella pertussis</i>	Vaccines	A											●																										
<i>Candida albicans</i>	Various												●					●																					
<i>Clostridium sp.</i>	Various												●	●	●	●	●																						
<i>Clostridium tetani</i>	Vaccines												●	●	●	●																							
<i>Corynebacterium diphtheriae</i>	Vaccines												●																										
<i>Escherichia coli</i>	rProtein	A	⊕ ⊕ ⊕ ⊕ ⊕ ⊕ ⊕					● ●		● *			●																										
<i>Haemophilus influenzae</i>	Vaccines												●		●																								
<i>Lactobacillus sp.</i>	Cultures												●	●			●																						
<i>Lactobacillus acidophilus</i>	Probiotics	A											●	●																									
<i>Lactobacillus plantarum</i>	Cultures	A											●	●	●	●	●																						
<i>Lactobacillus casei</i>	Probiotics	A											●	●	●	●	●																						
<i>Lactobacillus bulgaricus</i>	Cultures	A											●	●	●	●	●																						
<i>Lactococcus lactis sp.</i>	Various	A											●	●	●	●	●																						
<i>Leuconostoc mesenteroides</i>	Cultures	A						●					●	●	●	●	●			●	●	●	●	●	●	●													
<i>Lysinibacillus sphaericus</i>	Biopesticides												●	●																									
<i>Micromonospora sp.</i>	Antibiotics																																						
<i>Mycobacterium sp.</i>	Vaccines																																						
<i>Neisseria meningitidis</i>	Vaccines							●					●																										
<i>Nocardia sp.</i>	Antibiotics												●							●																			
<i>Pasteurella sp.</i>	Vaccines												●	●																									
<i>Pediococcus sp.</i>	Cultures												●	●	●	●	●			●	●	●	●	●	●	●													
<i>Pediococcus pentosaceus</i>	Cultures	A																																					
<i>Penicillium sp.</i>	Antibiotics																																						
<i>Pichia pastoris</i>	rProtein																																						
<i>Propionibacterium sp.</i>	Various							●	●				●							●	●	●	●																
<i>Pseudomonas sp.</i>	Various												●							●																			
<i>Saccharomyces cerevisiae</i>	rProtein		⊕ ⊕																																				
<i>Salmonella sp.</i>	Vaccines																			●																			
<i>Shigella sp.</i>	Vaccines																																						
<i>Staphylococcus sp.</i>	Vaccines																			●																			
<i>Staphylococcus carnosus</i>	Cultures	A											●	●	●	●	●			●	●	●	●	●	●														
<i>Streptococcus sp.</i>	Various												●							●																			
<i>Streptococcus thermophilus</i>	Cultures	A											●	●	●					●			●	●															
<i>Streptomyces sp.</i>	Various	A		⊕									●	●	●	●			●	●	●	●																	
<i>Trichoderma</i>	Enzymes																		●	●	●																		
<i>Vibrio cholerae</i>	Vaccines																			●																			
<i>Xanthomonas campestris</i>	Xanthan gum												●	●	●	●																							

Chart Key: ● Effective ⊕ Most effective * see product info sheet

Notes: ** Kosher and/or IPS version of this product is available

A - Application note available

For more specific strain recommendations, please contact your Account Director.

PRODUCT SELECTION: DIAGNOSTICS

		Hy-Peptone TSA/TSB ACF	Hy-Soy™	N-Z-Soy™ BL4	Hy-Peptone TSA/CASO	N-Z-Case™ Plus	N-Z-Amine™ AS	N-Z-Case™ M	N-Z-Case™ ME	Pepticas™	Edamini™ F	Peptonized Milk	Primateone™ RL
<i>B. subtilis</i> ³⁾	Count ⁴⁾	++	+++	+	+++	+++	+++	+++	+++	+++	+++	++	++
	Size ⁴⁾	+++	+++	+	++	+++	+	+	+	+	+	+	++
<i>E. coli</i>	Count	+++	+++	+++	+++	+++	+++	+++	++	+++	+++	++	+++
	Size	+++	++	+++	++	++	++	++	++	+	++	++	+++
<i>E. faecalis</i>	Count	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++
	Size	+++	+++	±	++	++	+	±	±	±	+	+	+
<i>C. albicans</i>	Count	+++	+++	++	+++	+++	+++	++	++	+++	+	+++	+++
	Size	+++	++++	+	+++	+++	++	+	+	±	++	++	++
<i>P. aeruginosa</i>	Count	+++	+++	+++	+++	+++	+++	+++	++	+++	+++	+++	++++
	Size	+++	+++	++++	+++	+++	++	++	++	++	+++	++	+++
<i>S. typhimurium</i>	Count	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++
	Size	+++	++	++	++	++	++	+	++	+	++	++	+++
<i>S. aureus</i>	Count	++	+	++	+++	+++	++	+++	+++	+++	++	+++	+++
	Size	+++	+	±	++	+++	+	+	+	±	+	++	+++
<i>S. epidermidis</i>	Count	+++	++	++	+++	+++	++	+	+	+	++++	+++	+++
	Size	+++	++	+	+++	++	+	±	±	±	++	++	+++
<i>S. pyogenes</i>	Count	+++	+++	++	+++	+++	++++	++++	+++	+++	+++	+++	++
	Size	++	++	+	++	++	+	+	+	+	++	++	+
<i>S. agalactiae</i>	Count	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++	+++
	Size	+++	+++	+	++	++	+	+	+	+	++	++	+
<i>S. pneumoniae</i>	Count	++++	+++	+++	+++	+++	+++	+++	+++	+++	++++	+++	+++
	Size	+++	+++	++	++	+++	++	++	+++	++	+++	++	++

CASO-agar (Casein-Soy Agar a.k.a. TSA; Tryptic Soy Agar) was used as reference medium; per test the typical colony diameter and total plate count of the CASO plates were set to 100% for the separate test organisms. The scores are based on a minimum of 8 separate measurements. Samples were tested at a concentration of 2% (with 0.5% NaCl 1.25% Agar Agar; pH adjusted to 7.3 before autoclaving for 15 minutes at 121°C). 1) *Bacillus subtilis* ATCC6633, *Escherichia coli* ATCC25922, *Enterococcus faecalis* ATCC29212, *Pseudomonas aeruginosa* ATCC27853, *Salmonella typhimurium* ATCC14028 and *Staphylococcus aureus* ATCC25923 were cultivated aerobically at 37°C for 24 hours. *Candida albicans* ATCC10231 and *Staphylococcus epidermidis* ATCC12228 were cultivated aerobically at 37°C for 48 hours. *Streptococcus pyogenes* ATCC19615, *Streptococcus pneumoniae* ATCC6305 and *Streptococcus agalactiae* ATCC12386 were cultivated under CO₂ atmosphere at 37°C for 48 hours on agar plates supplemented with 5% defibrinated sheep blood). 2) - represents no detected growth, ± represents moderate growth (colony size or count up to 25% relative to the reference plates), + represents fair growth (25 to 50%), ++ represents good growth (50 to 80%), +++ represents very good growth (80 to 120%) and ++++ represents excellent growth (more than 120%).

Bio-science product selection *continued*

NON-ANIMAL CHEMICAL AND PHYSICAL PROPERTIES: TYPICAL VALUES

Product	Source	Digestion	AN (%)	TN (%)	AN/TN ratio	Ash (%)	Moisture (%)	Molecular Weight Distribution (%)					
								>10 kDa	5-10 kDa	2-5 kDa	1-2 kDa	1k-500 Da	<500Da
Supplement Systems													
Sheff-CHO ACF	Various non animal	Enzymatic	---	11	---	8	4.1	26.1	6.5	7.7	10.4	12.5	36.8
Sheff-CHO Plus ACF	Various non animal	Enzymatic	---	7.5	---	6.9	3.5	16.6	5.2	7.4	12	14.8	44.1
Sheff-CHO PF ACF	Various non animal	Enzymatic	---	10.2	---	9.5	4	0	0.6	5.1	14.6	21.1	58.6
Sheff-CHO Plus PF ACF	Various non animal	Enzymatic	---	6.2	---	6.6	3	0	0.9	5.3	13.9	19.6	60.3
Sheff-Pulse I	Various non animal	Enzymatic	---	6.7	---	8.3	6	13.1	6.9	8.1	11.4	13.5	47.0
Sheff-Pulse II	Various non animal	Enzymatic	---	3.5	---	4.5	6	0	0	1.7	9.3	18.1	70.9
Sheff-Vax ACF	Various non animal	Enzymatic	---	9.8	---	10.5	3.4	27.9	5.6	6.3	10.2	13.1	36.9
Sheff-Vax Plus ACF	Various non animal	Enzymatic	---	8	---	6.4	3.2	43.2	6.2	6.4	8.3	9.5	26.4
Sheff-Vax PF ACF	Various non animal	Enzymatic	---	8.8	25	11.6	2.8	0.0	0.2	4.1	14.5	23.5	57.7
Sheff-Vax Plus PF ACF	Various non animal	Enzymatic	---	6.5	---	7.1	4.2	0	0	2.7	15.3	21.6	60.4
Sheff-Vax Plus PF ACF VP	Various non animal	Enzymatic	---	2.9	---	3.2	2	0	0.01	3.7	15.6	21.8	58.8
Hy-Express™ Systems													
Hy-Express™ System I	Various non animal	Enzymatic	3.0	12.0	24.2	9.0	3.7	0.4	2.0	10.4	17.3	20.8	49.1
Hy-Express™ System II	Various non animal	Enzymatic	3.0	13.0	25.0	5.0	4.0	0.0	0.1	2.5	10.4	18.2	68.9
Hy-Express™ System IV	Various non animal	Enzymatic	1.9	9.8	19.7	9.0	6.0	0.0	0.1	4.3	15.2	21.1	59.3
Hy-Express™ System V	Various non animal	Enzymatic	2.4	8.2	29.5	5.8	6.0	0.0	0.0	3.1	14.5	21.4	61.0
Hy-Express™ System VII	Various non animal	Enzymatic	4.9	11.4	43.2	9.2	5.5	0.0	0.8	6.5	14.7	17.0	61.0

Product	pH	Minerals (%)					
		Calcium	Iron	Magnesium	Phosphorus	Potassium	Sodium
Supplement Systems							
Sheff-CHO ACF	7.1	0.05	0.0040	0.257	0.46	1.9	2.4
Sheff-CHO Plus ACF	6.5	0.03	0.0026	0.186	0.47	1.67	1.63
Sheff-CHO PF ACF	7.1	0.058	0.0035	0.299	0.49	2.22	2.45
Sheff-CHO Plus PF ACF	6.8	0.04	0.0032	0.217	0.47	1.78	1.69
Sheff-Pulse I	6.2	0.07	0.0012	0.481	1.03	1.95	1.63
Sheff-Pulse II	5.6	0.052	0.0000	0.266	0.67	1.3	0.04
Sheff-Vax ACF	7.2	0.0289	0.0034	0.152	0.35	2.51	3.13
Sheff-Vax Plus ACF	7.2	0.018	0.0039	0.0914	0.24	1.45	2.06
Sheff-Vax PF ACF	7.7	0.02	0.0040	0.142	0.41	2.52	3.44
Sheff-Vax Plus PF ACF	7.3	0.0546	0.0033	0.311	0.52	2.22	2.51
Sheff-Vax Plus PF ACF VP	7.3	0.0123	0.00297	0.0398	0.123	0.868	1.13
Hy-Express™ Systems							
Hy-Express™ System I	4.9	0.04	0.009	0.2	1.04	2.2	1.20
Hy-Express™ System II	6.6	0.04	0.0018	0.061	0.951	0.172	1.68
Hy-Express™ System IV	6.2	0.10	0.0053	0.224	1.01	3.02	2.17
Hy-Express™ System V	6.1	0.04	0.0011	0.0329	0.154	0.173	0.946
Hy-Express™ System VII	6.3	0.03	0.005	0.0597	0.836	3.0	1.28

Bio-science product selection *continued*

NON-ANIMAL TOTAL AMINO ACIDS (MG/G): TYPICAL VALUES

Product	Alanine	Arginine	Aspartic Acid	Cysteine	Glutamic Acid	Glycine	Histidine	Isoleucine	Leucine	Lysine
Supplement Systems										
Sheff-CHO ACF	28	46	54	2	156	20	17	16	54	37
Sheff-CHO Plus ACF	17	28	36	1	106	14	10	13	29	21
Sheff-CHO PF ACF	24	53	64	1	191	25	15	23	44	30
Sheff-CHO Plus PF ACF	15	29	39	1	109	15	9	14	25	18
Sheff-Pulse I	18	37	39	N/A	88	12	9	11	24	21
Sheff-Pulse II	7	18	16	1	39	5	4	5	9	6
Sheff-Vax ACF	21	30	49	2	87	16	12	14	36	31
Sheff-Vax Plus ACF	26	30	57	1	100	15	13	13	42	43
Sheff-Vax PF ACF	25	37	66	3	112	24	13	22	38	32
Sheff-Vax Plus PF ACF	20	27	41	0	67	20	11	19	27	26
Sheff-Vax Plus PF ACF VP	8	14	25	0	45	10	5	8	15	13
Hy-Express™ Systems										
Hy-Express™ System I	43	53	84	0	176	27	15	26	46	58
Hy-Express™ System II	25	24	37	N/A	240	20	11	23	44	22
Hy-Express System IV	24	27	59	0	104	16	9	15	28	32
Hy-Express System V	26	38	64	0.7	85	16	16	24	41	30
Hy-Express System VII	24	28	43	0.7	73	19	10	21	35	26

Product	Methionine	Phenylalanine	Proline	Serine	Threonine	Tyrosine	Valine	Total Amino Acids	Total Free Amino Acids
Supplement Systems									
Sheff-CHO ACF	26	36	37	28	23	25	80	685	50
Sheff-CHO Plus ACF	6	22	28	21	15	14	42	423	59
Sheff-CHO PF ACF	10	32	46	36	23	22	73	712	58
Sheff-CHO Plus PF ACF	5	18	25	21	14	12	49	418	65
Sheff-Pulse I	4	19	16	18	13	9	58	396	92
Sheff-Pulse II	2	8	8	9	5	4	30	176	63
Sheff-Vax ACF	6	19	21	23	19	17	50	453	67
Sheff-Vax Plus ACF	5	26	21	24	20	19	37	492	42
Sheff-Vax PF ACF	7	24	27	29	22	20	24	525	68
Sheff-Vax Plus PF ACF	13	22	23	24	20	19	126	505	46
Sheff-Vax Plus PF ACF VP	3	9	10	13	8	1	79	266	25
Hy-Express™ Systems									
Hy-Express™ System I	9	28	34	37	28	22	92	778	41
Hy-Express™ System II	9	33	68	32	20	20	136	763	72
Hy-Express™ System IV	59	18	19	24	18	15	15	482	145
Hy-Express™ System V	11	30	25	26	21	22	51	528	52
Hy-Express™ System VII	8	24	18	23	20	30	65	468	229

Bio-science product selection *continued*

NON-ANIMAL CHEMICAL AND PHYSICAL PROPERTIES: TYPICAL VALUES

Product	Source	Digestion	AN (%)	TN (%)	AN/TN Ratio	Ash (%)	Moisture (%)	Molecular Weight Distribution (%)				
Ultra-filtered, Non-Animal Hydrolysates												
HyPep™ 1510	Soy	Enzymatic	2.0	9.2	21.7	9.7	3.6	0.0	0.2	4.6	14.3	22.8
HyPep™ 1511	Soy	Enzymatic	2.7	13.6	20.0	8.6	2.6	0.1	1.4	10.1	19.7	25.0
HyPep™ 1512	Soy	Enzymatic	2.2	8.8	25	11.6	2.8	0.0	0.1	3.2	13.5	21.1
HyPep™ 4601N	Wheat	Enzymatic	2.7	13.4	20.1	2.4	4.2	.17	1.25	5.5	13.1	21.1
HyPep™ 5603N	Rice	Enzymatic	2.3	12.8	18.0	3.9	3.8	0.5	1.6	9.7	19.7	24.6
HyPep™ 7504	Cotton	Enzymatic	1.9	9.9	19.5	10.7	3.2	0.1	2.0	9.9	19.1	21.9
UltraPep™ Soy	Soy	Enzymatic	2.5	8.0	31.3	11.0	5.2	0.2	1.4	5.1	12.0	18.4
Ultra-Filtered Yeast Extracts												
HyPep™ YE	Yeast	Enzymatic	5	10	50	9	4.3	0.0	0.1	3.4	10.2	13.7
UltraPep™ YE	Yeast	Enzymatic	4.4	9.6	45.8	13	2.9	0.0	0.6	6.0	14.7	19.3
Hy-Peptone												
Hy-Peptone TSA/TSB (ACF)	Various non animal	Enzymatic	3.0	10.0	30.0	15.0	6.0	0.2	1.2	6.9	15.4	19.4
Hy-Peptone 1009-K	Various	Enzymatic	1.6	5.6	28.1	8.7	6.0	9.3	10.5	12.4	12.2	11.1
Non-Animal Hydrolysates												
Amisoy™	Soy	Acid	8.6	13.0	66.0	10.1	3.7	0.0	0.0	0.4	4.5	18.5
Amisoy™ BF	Soy	Acid	6.7	9.3	72.1	14.0	5.0	0.0	0.0	0.9	7.3	16.5
Hy-Soy™	Soy	Enzymatic	1.9	9.1	21.0	10.7	2.5	0.0	0.3	5.3	16.4	23.5
Hy-Soy™ T	Soy	Enzymatic	1.8	7.9	22.8	6.6	3.9	0.0	0.2	5.0	16.3	23.8
N-Z-Soy™ BL4	Soy	Enzymatic	2.7	14.3	18.9	8.2	3.3	0.5	3.1	13.0	20.5	25.2
Hy-Express™ Peptone I	Various	Enzymatic	2.0	14.0	15.9	5.3	3.0	0.2	1.8	8.1	15.1	21.0
Hy-Express™ Peptone II	Various	Enzymatic	2.0	12.0	16.8	5.0	6.0	0.1	1.4	8.7	17.4	21.1
Hy-Pea™ 7404	Pea	Enzymatic	3.0	13.2	23.0	6.1	4.0	0.1	2.6	16.8	27.1	23.6
Hy-Pea™ 7477	Pea	Enzymatic	1.1	12.4	9.0	6.2	6.8	7.5	14.7	25.3	20.5	11.9
Hy-Rice™ 5377	Rice	Enzymatic	1.7	12.8	13.5	12.4	5.9	0.0	2.3	15.4	22.8	20.1
Yeast Products												
Hy-Yest™ 101	Yeast	Enzymatic	1.1	8.0	13.9	6.7	3.3	4.0	6.5	9.0	7.9	10.1
Hy-Yest™ 412	Yeast	Enzymatic	4.9	10.9	44.6	13.0	2.7	0.2	0.1	3.6	14.4	17.2
Hy-Yest™ 413	Yeast	Enzymatic	3.6	8.7	41.3	8.8	2.8	0.0	0.0	3.2	15.1	18.4
Hy-Yest™ 444	Yeast	Enzymatic	5.5	11.3	48.7	15.3	4.1	0.0	0.1	3.5	13.1	16.3
Hy-Yest™ 455	Yeast	Enzymatic	5.2	11.6	44.9	15.3	4.2	0.0	0.1	5.8	15.9	15.1
Hy-Yest™ 466	Yeast	Enzymatic	5	11.9	42.4	14.1	4.2	0.0	0.1	5.5	15.8	15.1
Hy-Yest™ 502	Yeast	Enzymatic	2.7	8.9	30.0	15.9	3.1	0.0	0.7	6.3	13.8	18.1
Hy-Yest™ 503	Yeast	Enzymatic	3.0	9.7	31.4	11.0	2.9	0.3	0.5	6.2	14.3	19.1
Hy-Yest™ 504	Yeast	Enzymatic	4.1	10.6	39.1	11.0	3.8	0.3	0.6	6.2	14.0	18.0
Hy-Yest™ 555	Yeast	Enzymatic	3.5	10.8	32.7	14.8	4.4	0.3	1.3	10.5	18.7	17.2
Sheftone™ Series												
Sheftone™ R	Various	Minimal	1.1	8.8	12.8	8.0	5.4	40.1	10.7	8.3	7.2	5.9

Product	pH	Minerals (%)					
		Calcium	Iron	Magnesium	Phosphorus	Potassium	Sodium
Ultra-filtered, Non-Animal Hydrolysates							
HyPep™ 1510	7.2	0.14	0.0004	0.28	0.03	3.40	3.30
HyPep™ 1511	7.0	0.03	0.0	0.02	0.25	0.70	4.04
HyPep™ 1512	7.7	0.02	0.004	0.142	0.41	2.52	3.44
HyPep™ 4601N	6.5	0.04	0.0014	0.039	0.17	0.16	1.03
HyPep™ 5603N	6.3	0.13	0.0002	0.003	0.21	2.10	1.20
HyPep™ 7504	7.1	0.03	0.004	0.15	2.86	3.47	3.79
UltraPep™ Soy	6.4	0.06	0.003	0.31	1.05	3.63	1.68
Ultra-Filtered Yeast Extracts							
HyPep™ YE	6.1	0.02	0.0061	0.051	1.02	3.21	1.50
UltraPep™ YE	6.7	0.25	0.0016	0.27	1.63	3.56	2.33
Hy-Peptone							
Hy-Peptone TSA/TSB (ACF)	7.0	0.04	0.0028	0.141	0.93	3.2	2.8
Hy-Peptone 1009-K	4.4	0.12	0.0107	0.59	1.5	2.61	0.35
Non-Animal Hydrolysates							
Amisoy™	5.8	2.00	0.0009	0.16	0.01	0.14	0.20
Amisoy™ BF	5.9	2.00	-	0.16	-	0.14	0.20
Hy-Soy™	7.1	0.20	0.006	0.25	0.27	3.00	2.70
Hy-Soy™ T	6.0	0.37	0.008	0.29	0.75	2.36	1.78
N-Z-Soy™ BL4	4.8	1.37	0.003	0.04	0.14	0.15	0.93
Hy-Express™ Peptone I	5.1	0.30	0.00018	0.004	0.023	0.04	0.13
Hy-Express™ Peptone II	5.2	0.06	0.0094	0.048	0.51	0.17	1.58
Hy-Pea™ 7404	4.7	0.30	0.02	0.02	0.69	0.40	1.70
Hy-Pea™ 7477	7.1	0.14	0.017	0.050	1.26	0.15	1.55
Hy-Rice™ 5377	6.3	0.06004	0.00215	0.03122	0.892	0.6846	1.227
Yeast Products							
Hy-Yest™ 101	6.4	0.15	0.0262	0.24	1.31	2.72	0.04
Hy-Yest™ 412	5.6	0.20	0.0004	0.03	3.10	3.60	1.70
Hy-Yest™ 413	7.1	0.12	0.0124	0.16	1.98	2.34	0.89
Hy-Yest™ 444	6.7	0.02	0.003	0.13	3.23	3.10	3.58
Hy-Yest™ 455	5.5	0.21	0.074	0.20	1.25	2.21	2.44
Hy-Yest™ 466	5.6	0.15	0.064	0.15	0.96	1.82	1.51
Hy-Yest™ 502	5.3	0.22	0.0156	0.268	1.89	3.17	4.59
Hy-Yest™ 503	5.1	0.17	0.0146	0.275	1.74	3.47	1.01
Hy-Yest™ 504	5.2	0.144	0.01	0.262	2.15	3.19	1.03
Hy-Yest™ 555	5.6	0.25	0.052	0.16	1.00	2.20	1.83
Sheftone™ Series							
Sheftone™ R	4.9	N/A	N/A	N/A	N/A	N/A	N/A

Bio-science product selection *continued*

NON-ANIMAL TOTAL AMINO ACIDS (MG/G): TYPICAL VALUES

Product	Alanine	Arginine	Aspartic Acid	Cysteine	Glutamic Acid	Glycine	Histidine	Isoleucine	Leucine	Lysine
Ultra-filtered, Non-Animal Hydrolysates										
HyPep™ 1510	26	41	73	4	122	25	14	23	40	38
HyPep™ 1511	30	62	92	1	166	29	17	28	53	52
HyPep™ 1512	25	37	66	3	112	24	13	22	38	32
HyPep™ 4601N	21	27	24	11	319	27	15	32	57	11
HyPep™ 5603N	48	72	81	8	174	37	20	34	66	35
HyPep™ 7504	27	71	58	4	113	26	14	18	33	21
UltraPep™ Soy	19	30	51	0	79	13	9	19	30	33
Ultra-Filtered Yeast Extracts										
HyPep™ YE	40	30	69	0	113	23	10	30	44	54
UltraPep™ YE	35	42	64	0	127	22	10	24	40	44
Hy-Peptone										
Hy-Peptone TSA/TSB (ACF)	35	37	62	1	186	23	11	24	41	39
Hy-Peptone 1009-K	20	14	21	0	48	12	6	10	22	11
Non-Animal Hydrolysates										
Amisoy™	51	66	82	N/A	176	33	15	36	61	54
Amisoy™ BF	38	46	69	N/A	126	23	12	30	49	45
Hy-Soy™	23	37	74	1	123	22	13	20	38	39
Hy-Soy™ T	25	44	74	1	119	20	13	21	41	42
N-Z-Soy™ BL4	30	65	90	0	170	35	17	25	50	106
Hy-Express™ Peptone I	27	50	65	N/A	259	25	17	29	56	42
Hy-Express™ Peptone II	23	48	55	0	221	25	15	24	51	34
Hy-Pea™ 7404	38	69	82	4	140	28	17	24	51	59
Hy-Pea™ 7477	27	61	80	0	134	25	15	32	56	53
Hy-Rice™ 5377	42	65	72	0	148	40	20	27	54	32
Yeast Products										
Hy-Yest™ 101	35	23	45	2	81	20	9	39	31	32
Hy-Yest™ 412	49	31	63	3	102	27	12	29	43	46
Hy-Yest™ 413	28	24	46	0	71	15	8	20	30	40
Hy-Yest™ 444	45	29	64	0	138	27	12	29	42	44
Hy-Yest™ 455	40	21	45	0	72	19	9	23	33	37
Hy-Yest™ 466	47	23	52	0	81	23	11	26	39	43
Hy-Yest™ 502	40	24	54	3	100	25	10	20	32	37
Hy-Yest™ 503	41	27	59	5	114	27	10	22	34	41
Hy-Yest™ 504	41	27	56	1	104	21	10	22	33	45
Hy-Yest™ 555	50	24	48	0	119	23	8	18	30	38

Product	Methionine	Phenylalanine	Proline	Serine	Threonine	Tyrosine	Valine	Total Amino Acids	Total Free Amino Acids
Ultra-filtered, Non-Animal Hydrolysates									
HyPep™ 1510	7	26	28	31	24	21	25	568	45
HyPep™ 1511	10	33	37	41	29	26	65	767	28
HyPep™ 1512	7	24	27	29	22	20	24	525	68
HyPep™ 4601N	12	46	115	38	21	29	32	835	63
HyPep™ 5603N	15	43	38	44	31	43	49	838	14
HyPep™ 7504	7	30	20	26	20	19	26	533	34
UltraPep™ Soy	2	20	39	21	17	17	19	418	96
Ultra-Filtered Yeast Extracts									
HyPep™ YE	9	26	24	30	28	14	99	641	319
UltraPep™ YE	8	25	21	29	24	16	72	601	208
Hy-Peptone									
Hy-Peptone TSA/TSB (ACF)	8	26	33	31	23	18	53	651	177
Hy-Peptone 1009-K	4	11	19	13	10	8	133	363	77
Non-Animal Hydrolysates									
Amisoy™	10	42	45	22	18	30	88	828	486
Amisoy™ BF	8	29	31	18	16	19	56	615	424
Hy-Soy™	5	26	29	32	22	18	22	542	63
Hy-Soy™ T	8	27	28	29	22	21	76	609	35
N-Z-Soy™ BL4	9	34	39	46	27	25	75	843	20
Hy-Express™ Peptone I	10	40	76	42	26	28	130	919	19
Hy-Express™ Peptone II	8	41	61	38	22	27	69	762	30
Hy-Pea™ 7404	6	32	28	40	28	26	27	699	7
Hy-Pea™ 7477	6	37	28	38	23	25	34	676	4
Hy-Rice™ 5377	15	37	34	53	32	40	39	751	12
Yeast Products									
Hy-Yest™ 101	6	18	17	23	23	14	23	441	54
Hy-Yest™ 412	9	25	22	28	28	21	35	573	255
Hy-Yest™ 413	5	18	36	19	19	18	23	421	86
Hy-Yest™ 444	9	24	23	26	27	18	33	590	299
Hy-Yest™ 455	9	20	16	22	23	15	25	429	200
Hy-Yest™ 466	10	23	18	26	26	18	29	495	295
Hy-Yest™ 502	7	18	19	20	21	15	25	470	107
Hy-Yest™ 503	5	20	20	26	27	14	26	518	95
Hy-Yest™ 504	7	22	21	24	23	35	27	521	142
Hy-Yest™ 555	6.5	17	16	24	24	12	22	479	133

Bio-science product selection *continued*

DAIRY & ANIMAL CHEMICAL AND PHYSICAL PROPERTIES: TYPICAL VALUES

Product	Source	Digestion	AN (%)	TN (%)	AN/TN Ratio	Ash (%)	Moisture (%)	Molecular Weight Distribution (%)					
								>10 kDa	5-10 kDa	2-5 kDa	1-2 kDa	1k-500 Da	<500Da
Supplement Systems													
Sheff-Vax PB10	Various	Enzymatic	3	7	---	32	6	0	0.5	3.9	9.3	17.2	69.1
Hy-Peptone													
Hy-Peptone TSA-CASO	Various	Enzymatic	4.4	12.0	36.8	9.8	6.0	0.0	0.4	3.6	11.6	23.7	60.7
Dairy Hydrolysates													
Amicase™	Casein	Acid	9.5	13.4	70.5	1.9	2.8	0.0	0.0	0.1	3.8	22.9	73.3
Hy-Case™ Amino	Casein	Acid	6.6	8.5	77.0	34.9	2.4	0.0	0.0	0.0	3.0	20.8	76.1
Hy-Case™ M	Casein	Acid	5.5	8.0	68.8	37.3	2.4	0.0	0.0	0.1	5.5	21.3	73.1
Hy-Case™ SF	Casein	Acid	10.6	13.6	77.6	0.4	2.9	0.0	0.0	0.0	2.7	30.2	67.1
N-Z-Amine™ A	Casein	Enzymatic	6.4	13.1	49.0	5.4	3.4	0.0	0.0	0.7	5.3	17.3	76.7
N-Z-Amine™ AS	Casein	Enzymatic	6.2	13.4	46.5	5.2	3.5	0.0	0.0	0.7	6.2	20.6	72.4
N-Z-Amine™ EKC	Casein	Enzymatic	3.7	13.3	28.0	5.1	4.5	0.0	0.0	1.0	6.6	20.5	71.9
N-Z-Case™ M	Casein	Enzymatic	4.9	13.4	36.2	6.5	3.1	0.0	0.0	2.4	9.8	28.6	59.2
N-Z-Case™ ME	Casein	Enzymatic	4.4	13.8	32.2	5.5	3.4	0.0	0.2	6.4	12.5	26.5	54.5
N-Z-Case™ Plus	Casein	Enzymatic	6.5	13.3	49.3	5.6	3.3	0.0	0.1	4.0	7.6	17.8	70.5
N-Z-Case™ TT	Casein	Enzymatic	5.3	13.3	39.6	5.6	3.2	0.0	0.0	1.5	7.6	24.1	66.9
Pepticase™	Casein	Enzymatic	4.1	13.2	31.4	9.7	3.3	0.0	0.0	4.2	12.0	27.5	56.3
Edamin™ F	Lactalbumin	Enzymatic	4.7	12.6	37.7	5.5	2.8	0.0	0.0	2.3	11.1	25.4	61.2
Peptonized Milk	Milk Solids	Enzymatic	1.7	5.7	29.8	7.5	3.9	0.0	0.2	4.6	11.7	23.8	59.7
Animal Hydrolysates													
Primatone™ HS	Meat	Enzymatic	4.6	12.5	36.8	12.3	2.6	2.2	4.0	9.7	11.4	15.1	57.6
Primatone™ HS/UF	Meat	Enzymatic	5.1	12.2	41.9	12.8	3.3	0.2	1.4	6.9	11.8	16.7	63.1
Primatone™ RL	Meat	Enzymatic	5.6	12.1	46.4	13.5	2.8	0.3	1.1	3.4	6.5	16.0	72.8
Primatone™ RL/UF	Meat	Enzymatic	5.9	11.8	49.6	13.8	2.6	0.0	0.3	2.2	5.8	15.7	76.0
Sheftone™ Series													
Sheftone™ EK	Various	Minimal	1.4	5.9	23.1	10.0	4.0	11.0	11.0	9.0	10.6	9.8	48.5

Product	pH	Minerals (%)					
		Calcium	Iron	Magnesium	Phosphorus	Potassium	Sodium
Supplement Systems							
Sheff-Vax PB10	7.0	0.0252	0.0013	0.0198	1.52	3.0	9.32
Hy-Peptone							
Hy-Peptone TSA-CASO	6.9	0.03	0.0015	0.08	1.03	0.85	3.19
Dairy Hydrolysates							
Amicase™	6.6	0.03	-	0.002	0.03	0.01	0.99
Hy-Case™ Amino	5.7	0.03	0.001	0.03	-	0.02	14.40
Hy-Case™ M	6.9	0.02	0.001	0.01	0.48	0.07	13.70
Hy-Case™ SF	6.4	0.03	<0.0004	0.01	0.05	0.01	0.31
N-Z-Amine™ A	6.8	0.02	0.003	0.01	0.89	0.08	2.56
N-Z-Amine™ AS	6.8	0.04	0.004	0.02	0.97	0.08	2.49
N-Z-Amine™ EKC	6.8	0.04	0.002	0.02	0.76	0.22	1.69
N-Z-Case™ M	7.2	0.03	0.002	0.02	0.96	0.09	2.20
N-Z-Case™ ME	7.2	0.03	0.001	0.005	0.83	0.04	2.82
N-Z-Case™ Plus	6.8	0.03	0.002	0.1	0.82	0.13	1.40
N-Z-Case™ TT	7.2	0.02	0.003	0.006	0.76	0.09	2.30
Pepticase™	7.2	0.03	0.001	0.02	1.52	0.03	3.34
Edamin™ F	6.9	0.14	0.0029	0.035	0.52	0.41	2.26
Peptonized Milk	6.3	1.23	0.004	0.19	1.02	1.71	0.53
Animal Hydrolysates							
Primatone™ HS	7.3	0.03	0.004	0.04	0.65	1.50	4.92
Primatone™ HS/UF	7.1	0.03	0.002	0.03	0.69	1.39	5.15
Primatone™ RL	7.3	0.03	0.004	0.04	0.85	1.78	5.31
Primatone™ RL/UF	7.3	0.01	0.022	0.02	1.09	1.45	5.45
Sheftone™ Series							
Sheftone™ EK	4.9	0.46	0.0116	0.604	1.89	2.95	0.42

Bio-science product selection *continued*

DAIRY & ANIMAL TOTAL AMINO ACIDS (MG/G): TYPICAL VALUES

Product	Alanine	Arginine	Aspartic Acid	Cysteine	Glutamic Acid	Glycine	Histidine	Isoleucine	Leucine	Lysine
Supplement Systems										
Sheff-Vax PB10	32	31	51	3	119	32	16	30	55	51
Hy-Peptone										
Hy-Peptone TSA-CASO	23	36	59	0	167	20	21	36	62	56
Dairy Hydrolysates										
Amicase™	63	42	35	NA	175	22	27	57	88	104
Hy-Case™ Amino	25	22	44	NA	140	11	14	30	53	58
Hy-Case™ M	16	21	44	NA	135	10	14	27	49	54
Hy-Case™ SF	48	41	52	NA	151	18	25	52	92	109
N-Z-Amine™ A	23	31	58	0	170	13	19	40	76	72
N-Z-Amine™ AS	25	34	63	0	187	14	21	41	68	80
N-Z-Amine™ EKC	24	32	60	0	176	13	20	36	69	73
N-Z-Case™ M	24	35	59	1	180	15	21	41	71	73
N-Z-Case™ ME	26	35	67	1	196	14	21	42	69	81
N-Z-Case™ Plus	25	34	58	1	171	15	20	40	70	71
N-Z-Case™ TT	23	33	59	1	178	13	21	38	73	77
Pepticase™	24	35	60	1	178	15	18	38	73	74
Edamin™ F	26	24	68	0	139	12	14	38	66	68
Peptonized Milk	18	21	41	0	103	16	11	21	42	37
Animal Hydrolysates										
Primatone™ HS	51	45	73	1	117	57	17	26	59	66
Primatone™ HS/UF	43	29	60	0	106	45	14	23	53	111
Primatone™ RL	47	43	68	0	92	35	23	19	66	68
Primatone™ RL/UF	49	39	73	0	105	34	22	21	63	74

Product	Methionine	Phenylalanine	Proline	Serine	Threonine	Tyrosine	Valine	Total Amino Acids	Total Free Amino Acids
Supplement Systems									
Sheff-Vax PB10	16	27	56	32	27	14	38	630	215
Hy-Peptone									
Hy-Peptone TSA-CASO	19	40	71	43	32	21	94	800	210
Dairy Hydrolysates									
Amicase™	29	46	130	21	27	37	73	977	655
Hy-Case™ Amino	15	28	63	24	20	10	39	596	418
Hy-Case™ M	14	26	59	32	23	16	36	577	540
Hy-Case™ SF	26	42	117	36	33	20	87	948	703
N-Z-Amine™ A	17	37	122	40	32	49	48	847	414
N-Z-Amine™ AS	12	40	95	44	35	50	53	861	401
N-Z-Amine™ EKC	7	41	94	43	33	42	48	811	109
N-Z-Case™ M	17	40	88	45	35	17	58	819	234
N-Z-Case™ ME	10	43	95	47	36	19	51	854	186
N-Z-Case™ Plus	13	40	91	45	34	46	49	823	382
N-Z-Case™ TT	17	43	96	40	34	62	45	852	244
Pepticase™	12	42	93	44	33	87	47	873	132
Edamin™ F	16	27	55	37	39	17	145	791	256
Peptonized Milk	10	22	36	28	19	23	115	563	38
Animal Hydrolysates									
Primatone™ HS	13	28	44	34	30	14	24	698	261
Primatone™ HS/UF	12	27	35	31	27	16	110	743	280
Primatone™ RL	12	33	30	33	31	18	44	661	401
Primatone™ RL/UF	12	32	30	35	32	21	44	686	366

Appendix

Methods used to determine typical data

Degree of Hydrolysis

The degree of hydrolysis was estimated by the AN/TN ratio. Total Nitrogen (TN) values were determined using a Dumas or Kjehldahl method (2). The concentrations of alpha-amino nitrogen (AN) was determined by the formaldehyde titration method described in the U.S. Pharmacopoeia (1). The ratio of AN/TN is an estimate of the degree of hydrolysis of the protein.

Amino Acids

Total amino acids were determined with an HPLC after acid hydrolysis. With this method tryptophan is completely destroyed. Cysteine is partly destroyed and is reported as the dimer cystine. Glutamine and asparagine are converted into glutamic acid and aspartic acid respectively.

Ash

Ash values were measured after heating at 650 °C.

- > Details of each analytical method are available upon request.
- > The product specification does not necessarily include all above described analysis.

Moisture

The moisture level was determined as the loss on drying upon heating at 100-103 °C.

Minerals

Minerals were determined with ICP.

pH

The pH was measured in a 2% solution at 25 °C.

Molecular Weight Distribution

The molecular weight distribution was determined by HPLC/Size Exclusion Chromatography using acetonitrile/water eluent and absorbance at 220 nm for detection. It should be noted that other methods will yield different results therefore these values should only be used to compare products utilizing the same method. Molecular weight distribution was determined of the soluble fraction of the sample only.

List of microorganisms

Organism	Product	Application	References
<i>Absidia glauca</i>	Edamin™	oxygenation of steroids	15
<i>Achyla ambisexualis</i>	Edamin™	growth medium	123
<i>Achyla heterosexualis</i>	Edamin™	growth medium	119
<i>Actinomadura ferruginea</i>	N-Z-Amine™ A N-Z-Case™	production of antitumor antibiotic production of antitumor antibiotic	629 629
<i>Actinomadura hibisca</i>	N-Z-Amine™ A N-Z-Case™	production of antibiotics production of antibiotics	478, 483-484 478, 483-484
<i>Actinomadura kijaniata</i>	N-Z-Amine™ A	synthesis of antibiotics	223
<i>Actinomadura luzonensis</i>	N-Z-Amine™ A	antibiotic production	205
<i>Actinomadura madureae</i>	N-Z-Amine™ A	production of antibiotics	450, 476
<i>Actinomadura mellaura</i>	N-Z-Amine™	antibiotic production	357
<i>Actinomadura oligospora</i>	N-Z-Amine™ A	antibiotic production	434, 436
<i>Actinomadura sp</i>	N-Z-Amine™ N-Z-Amine™ A N-Z-Case™ Peptonized Milk	production of antibiotics antibiotic production antibiotic production antibiotic production	323 448, 482, 503 367 367
<i>Actinomadura verrucospora</i>	N-Z-Amine™ A N-Z-Case™	antibiotic production antibiotic production	366, 412 366, 412
<i>Actinomyces griseus</i>	N-Z-Amine™	production of Streptomycin	11, 13
<i>Actinomycete spp.</i>	N-Z-Amine™ A	growth medium	394
<i>Actinoplanacetace sp.</i>	Hy-Case™ SF	fermentation medium	507, 529, 558
<i>Actinoplanes sp.</i>	Hy-Case™ SF N-Z-Amine™ A	fermentation medium growth medium	540, 556 428
<i>Actinopolyspora sp.</i>	Amicase™ Hy-Soy™ N-Z-Amine™ A	growth medium, antibiotic production growth medium, antibiotic production growth medium, antibiotic production	620 620 620
<i>Actinosynnema mirum</i>	N-Z-Amine™	growth medium	175
<i>Aerobacter aerogenes</i>	N-Z-Amine™	hydrolysis of penicillins	63
<i>Aeromonas punctata</i>	Edamin™ S Hy-Case™ SF N-Z-Amine™ A N-Z-Amine™ AS Primatone™	production of L-serine dehydratase production of L-serine dehydratase production of L-serine dehydratase production of L-serine dehydratase production of L-serine dehydratase	154 154 154 154 154
<i>Agrobacterium tumefaciens</i>	N-Z-Case™	cobalamin production	463
<i>Alcaligenes faecalis</i>	N-Z-Amine™ N-Z-Amine™	hydrolysis of penicillins production of vitamin B12	63 30
<i>Alcaligenes sp</i>	N-Z-Amine™ A	kalimantacin production	589
<i>Amorphosporangium globisporum</i>	N-Z-Amine™ A	growth medium	272
<i>Amycolata sp.</i>	Amicase™ Hy-Soy™ N-Z-Amine™ A	growth medium, antibiotic production growth medium, antibiotic production growth medium, antibiotic production	620 620 620
<i>Aspergillus fumigatus</i>	N-Z-Amine™ A	production of farnesyltransferase inhibitor	584
<i>Aspergillus nidulans</i>	N-Z-Amine™ A N-Z-Case™	antibiotic production antibiotic production	240 240
<i>Aspergillus niger</i>	N-Z-Amine™ A	growth medium	392
<i>Aspergillus rugulosus</i>	N-Z-Amine™ A N-Z-Case™	production of antifungal agent production of antifungal agent	165 165
<i>Aspergillus spp.</i>	Edamin™ N-Z-Amine™	growth medium growth medium	25 25
<i>Aspergillus tamarii</i>	Hy-Soy™ T N-Z-Amine™ A	antibiotic production antibiotic production	147, 162 147, 162
<i>Aspergillus terreus</i>	N-Z-Amine™ Peptonized Milk	production of a hypocholesteremic production of a hypocholesteremic	225, 241, 275, 277 225, 241, 275, 277
<i>Aspergillus variecolor</i>	N-Z-Amine™ A	production of variecolin	527
<i>Aureobasidium pullulans</i>	N-Z-Amine™ A	growth medium	392
<i>Bacillus cereus</i>	Hy-Case™ SF	growth medium	183
<i>Bacillus pumilis</i>	N-Z-Case™ Pepticase™	production of antiviral agent production of peroxidase	480 530, 624
<i>Bacillus sp.</i>	Primatone™	production of alternanase	614
<i>Bacillus sphaericus</i>	Peptonized Milk	growth medium	412
<i>Bacillus subtilis</i>	Hy-Soy™ T N-Z-Amine™	growth medium growth medium	634, 655 309

List of microorganisms *continued*

Organism	Product	Application	References
<i>Bacillus thuringiensis</i>	Amisoy™ N-Z-Amine™ A	production of antibiotics production of antibiotics	550 313
<i>Bacteroides fragilis</i>	N-Z-Amine™ A	growth medium	398
<i>Bifidobacterium infantis</i>	N-Z-Amine™ A N-Z-Case™ Plus	growth medium growth medium	635 635
<i>Bifidobacterium longum</i>	N-Z-Amine™ A N-Z-Case™ Plus	growth medium growth medium	635 635
<i>Bortedella bronchiseptica</i>	N-Z-Amine™ AS	vaccine stabilization	383
<i>Bortedella pertussis</i>	Hy-Pep™ Hy-Soy™ N-Z-Amine™	production of toxin cultivating pathogenic bacteria hydrolysis of penicillins	689 618 63
<i>Botryosphaeria spp.</i>	Edamin™	hydroxylation of steroids	408
<i>Botyrodiploida theobromae</i>	Edamin™	hydroxylation of steroids	408
<i>Brevibacterium spp.</i>	N-Z-Amine™ Pepticase™	production of tryptophan mononitrilase production	413 395
<i>Brucella abortus</i>	N-Z-Amine™ A	growth medium	398
<i>Candida albicans</i>	N-Z-Amine™ A	growth medium	398
<i>Caryophanon latum</i>	Edamin™	growth medium	80
<i>Catenuloplanes japonicus</i>	N-Z-Amine™ A	growth medium	310
<i>Cephalosporium sp.</i>	Edamin™	production of cephalosporin	355
<i>Chainia spp.</i>	N-Z-Amine™ A	growth medium	428
<i>Chromobacterium violaceum</i>	N-Z-Amine™ A	antibiotic production	234, 257
<i>Citrobacter freundii</i>	N-Z-Amine™ A	synthesis of colicins	230
<i>Cladobotryum varium</i>	N-Z-Amine™ A	antibiotics production	662
<i>Clostridium acetobutylicum</i>	Hy-Case™	production of acetone and butanol	553
<i>Clostridium botulinum</i>	Edamin™ S Hy-Soy™ N-Z-Amine™ A	growth medium growth medium protease involved in toxin production	79 79 491
<i>Clostridium diphtheriae</i>	HyPep™ 5603	production of toxin	689
<i>Clostridium novyi</i>	N-Z-Amine™ A	growth medium	398
<i>Clostridium perfringens</i>	Edamin™ S Hy-Soy™ N-Z-Amine™ N-Z-Amine™ A N-Z-Amine™ YT N-Z-Case™	growth medium growth medium growth medium growth medium growth medium growth medium	79 79 401 398 43 393
<i>Clostridium sporogenes</i>	Edamin™ S Hy-Soy™	growth medium growth medium	79 79
<i>Clostridium tertium</i>	N-Z-Amine™ A	growth medium	38
<i>Clostridium tetani</i>	Amisoy™ HyPep™ 5603 HyPep™ 4601 Hy-Soy™ Hy-Soy™ T N-Z-Case™ N-Z-Case™ TT N-Z-Soy™ NZ-Soy N-Z-Soy™ BL 4 N-Z-Soy™ BL 7 Pepticase™	production of tetanus toxin in animal free media production of toxin production of toxin production of tetanus toxin in animal free media production of tetanus toxin in animal free media production of tetanus toxin production of tetanus toxin in animal free media production of tetanus toxin in animal free media cultivating pathogenic bacteria production of tetanus toxin in animal free media production of tetanus toxin in animal free media production of tetanus toxin	658 689 689 658, 687, 688, 690 658 7,12, 18, 27,49, 658 658, 687 658 618 658 658 658 7,12, 18, 27,49
<i>Conidiobolus coronatus</i>	Edamin™	synthesis of diphtheria toxins	552
<i>Coniothyrium fuckelii</i>	Amicase™	production of lovastin	577
<i>Corynebacterium diphtheriae</i>	Hy-Soy™ N-Z-Amine™ A	cultivating pathogenic bacteria synthesis of diphtheria toxins	618 1, 67, 111
<i>Corynebacterium fascians</i>	N-Z-Amine™	hydrolysis of penicillins	63
<i>Crithidia fasciculata</i>	N-Z-Amine™ AS	assay for biopterin	115, 124, 213
<i>Crithidia spp.</i>	N-Z-Amine™ AS	growth medium	117
<i>Dactylosporangium spp.</i>	N-Z-Amine™ A	growth medium	428
<i>Dermatophilus spp.</i>	N-Z-Amine™ A	growth medium	69
<i>Dictyostelium discoideum</i>	N-Z-Case™	growth medium	430
<i>Dictyostelium mucoroides</i>	N-Z-Case™	growth medium	430
<i>Diplococcus pneumoniae</i>	N-Z-Amine™ A	growth medium	398

<i>Elactomyces cylindrosporae</i>	N-Z-Case™	antibiotic production	369
<i>Empedobacter haloabium</i>	N-Z-Case™	antibiotic production	326
<i>Erwinia aroideae</i>	Hy-Soy™ T	synthesis of L-asparaginase	92
<i>Escherichia coli</i>	Amisoy™	production of glycerol	616
	Amisoy™	DNA vaccine production	648
	Edamin™ S	synthesis of galactosidase	390
	Hy-Case™ Amino	growth medium	390
	Hy-Case™ Amino	production of human growth hormone	504
	Hy-Case™ M	production of human growth hormone	504
	Hy-Case™ SF	production of lysine (arginine) decarboxylase	20-22, 50, 109
	Hy-Case™ SF	antibody production with plasmid host cells	676
	HyPep™ 4601	DNA vaccine production	648
	Hy-Rice™ 5303	DNA vaccine production	648
	Hy-Soy™	growth medium	429, 625-626
	Hy-Yest™ 412	growth medium for plasmid host cells	666,653
	N-Z-Amine™ A	production of swine growth hormone	452
	N-Z-Amine™ A	genetic studies	161, 177-178, 188-189, 203, 210, 229, 260-261, 265-270, 284, 288-290, 328, 372-374, 377, 416, 418, 426, 443, 445-446, 455-456, 460, 462, 464, 473, 489-490, 492-494, 496-497, 499
<i>Frankia spp.</i>	N-Z-Amine™ A	growth medium	596
<i>Fungal sp.</i>	N-Z-Amine™ A	synthesis of bovine growth hormone	449
	N-Z-Amine™ A	synthesis of colicins	230
	N-Z-Amine™ AS	growth medium in the production of polypeptides	651
	N-Z-Case™	genetic studies	66, 99
	N-Z-Case™	mutation studies	39
	N-Z-Case™ Plus	growth medium for plasmid host cells	666
	N-Z-Amine™ A	growth medium	254
	N-Z-Amine™ A	production of cellulases	414
	N-Z-Case™	production of cellulases	414
<i>Fusobacterium necrophorum</i>	N-Z-Amine™ A	growth medium	398
<i>Geodermatophilus obscurus</i>	N-Z-Amine™ A	growth medium	68
<i>Glycomyces harbinensis</i>	N-Z-Amine™ A	antibiotic production	409
<i>Haemophilus influenzae</i>	Hy-Soy™	growth medium	354
	N-Z-Amine™ A	growth medium	398
<i>Halobacterium halobium</i>	Hy-Case™ SF	growth medium	244
<i>Halobacterium saccharovorum</i>	Hy-Case™ SF	growth medium	244
<i>Halobacterium spp.</i>	Hy-Case™ SF	growth medium	96, 360, 399
<i>Haloferax denitrificans</i>	Hy-Case™ SF	Nitrite reductase preparation	597
<i>Helicobacter pylori</i>	Hy-Soy™	cultivating pathogenic bacteria	618
	N-Z-Soy™	cultivating pathogenic bacteria	618
<i>Helminthosporium sativum</i>	N-Z-Amine™ A	production of cochlorquinone	477
<i>Herpes virus (Marek's Disease)</i>	N-Z-Amine™ AS	stabilization during lyophilization	125
	N-Z-Amine™ AS	stabilization during lyophilization	125
<i>Herpetomonas spp.</i>	N-Z-Amine™ AS	growth medium	117
<i>Klebsiella pneumoniae</i>	Hy-Case™ SF	growth medium	361
<i>Lactobacillus acidophilus</i>	Peptonized Milk	growth medium	441
	Primatone™	growth medium	441
<i>Lactobacillus amylovorus</i>	Cottonseed	production of lactic acid	523
	Hy-Soy™	production of lactic acid	523
	N-Z-Soy™ BL	production of lactic acid	523
<i>Lactobacillus arabinosus</i>	Hy-Soy™	growth medium	43
	N-Z-Amine™ YT	growth medium	43
<i>Lactobacillus bifidus</i>	N-Z-Case™	growth medium	26
<i>Lactobacillus casei</i>	Edamin™ S	growth medium	652
<i>Lactobacillus fermentum</i>	N-Z-Case™	synthesis of glutamate dehydrogenase	348
<i>Lactobacillus lactis</i>	N-Z-Amine™	starter media for cheese cultures	468
	N-Z-Case™	synthesis of glutamate dehydrogenase	348
<i>Lactobacillus reuteri</i>	Hy-Soy™	growth medium in enzyme production	630
<i>Lasiodiplodia theobromae</i>	Edamin™	hydroxylation of steroids	408
<i>Leptomonas spp.</i>	N-Z-Amine™ AS	growth medium	117, 124
<i>Leuconostoc mesenteroides</i>	N-Z-Amine™ YT	growth medium	43
<i>Metarrhizium anisopliae</i>	N-Z-Amine™ A	production of enzymes	646
<i>Microbacterium sp.</i>	Hy-Soy™	culture medium	587
<i>Microbispora spp.</i>	N-Z-Amine™ A	growth medium	69

List of microorganisms *continued*

Organism	Product	Application	References
<i>Micromonospora carbonacea</i>	N-Z-Amine™ A	production of evernomicin	76
<i>Micromonospora chalcea</i>	N-Z-Amine™ A	antibiotic production	73, 134, 206
<i>Micromonospora citrea</i>	N-Z-Amine™ A	production of citreamicins	481
<i>Micromonospora coerulea</i>	N-Z-Amine™ A	growth medium	245
<i>Micromonospora danubiensis</i>	N-Z-Amine™ A	production of sisomicin	295-296
<i>Micromonospora echinospora</i>	N-Z-Amine™ A	antibiotic production	52, 56, 340
<i>Micromonospora fusca</i>	N-Z-Amine™ A	synthesis of gentamicins	52, 56
<i>Micromonospora griseorubida</i>	N-Z-Amine™ A	synthesis of mycinamicins	208
<i>Micromonospora halophytica</i>	N-Z-Amine™ A	production of halomicin	77
<i>Micromonospora heviziensis</i>	N-Z-Amine™ A	growth medium	311
<i>Micromonospora lacustris</i>	N-Z-Amine™ A	production of rifamycins	135
<i>Micromonospora megalomicea</i>	N-Z-Amine™ A	growth medium	102
<i>Micromonospora olivoasterospora</i>	N-Z-Amine™ A	growth medium	143
<i>Micromonospora polytrota</i>	N-Z-Amine™ A	growth medium	311
<i>Micromonospora purpurea</i>	N-Z-Amine™ A	production of gentamicins	144, 146, 166, 239, 300
<i>Micromonospora purpureochromogenes</i>	N-Z-Amine™ A	growth medium	420
<i>Micromonospora rhodorangea</i>	N-Z-Amine™ A	growth medium	83
<i>Micromonospora rosaria</i>	N-Z-Amine™ A	synthesis of roسامicin	118
	Peptidase™	production of dienone macrolides	97
		production of sisomicin	613
<i>Micromonospora rosea</i>	N-Z-Amine™	antibiotic production	279
<i>Micromonospora sagamiensis</i>	N-Z-Amine™ A	antibiotic production	170
<i>Micromonospora saitamica</i>	N-Z-Amine™ A	antibiotic production	167
<i>Micromonospora scalabitana</i>	N-Z-Amine™ A	synthesis of gentamicins	237
<i>Micromonospora viridifaciens</i>	N-Z-Amine™ A	antibiotic production	148
<i>Micropolyspora brevicatena</i>	N-Z-Amine™ A	growth medium	47
<i>Micropolyspora spp.</i>	N-Z-Amine™ A	growth medium	69
<i>Moraxella bovis</i>	N-Z-Amine™ A	vaccine production	358
<i>Mucor spp.</i>	Edamin™	oxygenation of steroids	15
<i>Mycobacterium avium</i>	N-Z-Case™	growth medium	36
<i>Mycobacterium spp.</i>	N-Z-Amine™	transformations of sterols	93
<i>Mycobacterium tuberculosis</i>	N-Z-Amine™ A	growth medium	69
	N-Z-Case™	growth medium	5, 36
		growth medium	36
<i>Mycobacterium vole</i>	N-Z-Case™	culture medium	53, 64
<i>Mycoplasma</i>	Hy-Soy™	growth medium	500
<i>Mycoplasma pneumoniae</i>	Hy-Soy™	vaccine stabilization	384
<i>Myxovirus influenza</i>	N-Z-Amine™ A	cultivating pathogenic bacteria	618
<i>Neisseria meningitidis</i>	Amisoy™	growth medium	398
	N-Z-Amine™ A	growth medium	246
<i>Nocardia asteroides</i>	N-Z-Amine™ A	growth medium	514
<i>Nocardia autotrophica</i>	Hy-Case™ SF	growth medium	246
<i>Nocardia brasiliensis</i>	N-Z-Amine™ A	growth medium	301
<i>Nocardia canicruria</i>	Hy-Case™	hydroxylation of steroids	301
	Hy-Soy™	hydroxylation of steroids	421
<i>Nocardia orientalis</i>	Amicase™	antibiotic production	421
	Hy-Case™	antibiotic production	421
	N-Z-Amine™ A	antibiotic production	421
<i>Nocardia sp.</i>	Hy-Yest™ 412	nocathiacin antibiotic production	645
	N-Z-Amine™ A	antibiotic production	228, 236, 479
	N-Z-Case™	manufacture of naphthofuranone compounds	640
<i>Nocardiopsis spp.</i>	N-Z-Amine™ A	growth medium	428
<i>Nodulisporium sp.</i>	Amicase™	growth medium	615
<i>Ochromonas danica</i>	N-Z-Amine™	growth medium	100
<i>Ochromonas malhamensis</i>	N-Z-Amine™	growth medium	100
<i>Oerskovia spp.</i>	N-Z-Amine™ A	growth medium	428
<i>Pasteurella haemolytica</i>	Hy-Case™ Amino	growth medium, vaccine production	382
	N-Z-Amine™ A	growth medium, vaccine production	382
<i>Pasteurella multocida</i>	Hy-Case™ Amino	growth medium, vaccine production	382
	N-Z-Amine™ A	growth medium, vaccine production	382

<i>Penicillium simplicissimum</i>	N-Z-Amine™ A	growth medium, antibiotic production	638
<i>Penicillium adametzioides</i>	Amicase™	production of compactin	602
<i>Penicillium levitum</i>	N-Z-Amine™ A	growth medium	392
<i>Penicillium sp.</i>	Edamin™	growth medium, hydroxylation of steroids	25
	N-Z-Amine™	growth medium, hydroxylation of steroids	25
	N-Z-Amine™ A	production of methionine-aminopeptidase inhibitor	637
	N-Z-Case™	growth medium	611
<i>Penicillium waksmanii</i>	N-Z-Amine™	growth medium	352
<i>Peniophora affinis</i>	N-Z-Amine™ A	antibiotic production	197
<i>Peptococcus aerogenes</i>	N-Z-Amine™ A	growth medium	398
<i>Peptostreptococcus sp.</i>	N-Z-Amine™ A	growth medium	398
<i>Phycomyces spp.</i>	Edamin™	oxygenation of steroids	15
<i>Pichia pastoris</i>	Hy-Yest™ 444	growth medium, production of recombinant proteins	675
<i>Polyangium brachysporum</i>	N-Z-Case™	antibiotic production	451
<i>Prevotella ruminicola</i>	Pepticase™	growth medium	639
<i>Promicromonospora spp.</i>	N-Z-Amine™ A	growth medium	428
<i>Proteus rettgeri</i>	N-Z-Amine™	hydrolysis of penicillins	63
<i>Proteus sphingidis</i>	N-Z-Amine™	hydrolysis of penicillins	63
<i>Proteus vulgaris</i>	Primatone™	synthesis of L-serine dehydratase	138
<i>Pseudomonas aeruginosa</i>	N-Z-Amine™ A	growth medium	398
<i>Pseudomonas cepacia</i>	N-Z-Amine™ A	antibiotic production	324, 410
<i>Pseudomonas delafieldii</i>	Edamin™	synthesis of L-serine dehydratase	173
	Hy-Soy™	synthesis of L-serine dehydratase	173
	Primatone™	synthesis of L-serine dehydratase	173
	Primatone™ HS	synthesis of L-serine dehydratase	173
<i>Pseudomonas denitrificans</i>	N-Z-Case™	cobalamin production	463
<i>Pseudomonas spp.</i>	N-Z-Amine™	production of vitamin B12	30
	N-Z-Amine™ A	production of antitumor substance	294
<i>Pseudomonas testosteroni</i>	Pepticase™	synthesis of glutaminase	129
<i>Pseudonocardia sp.</i>	Amicase™	growth medium, antibiotic production	620
	Hy-Soy™	growth medium, antibiotic production	620
	N-Z-Amine™ A	growth medium	428, 606
	N-Z-Amine™ A	growth medium, antibiotic production	620
<i>Pycnidiphora dispersa</i>	N-Z-Amine™ A	growth medium, antibiotic production	541
<i>Rhizopus oligosporus</i>	Hy-Soy™	growth medium	307
	N-Z-Amine™ A	growth medium	307
<i>Rhizopus spp.</i>	Edamin™	oxygenation of steroids	15
<i>Rhodococcus rhodochrous</i>	N-Z-Amine™	growth medium	498, 508
<i>Rhodococcus spp.</i>	N-Z-Amine™ A	growth medium	428
	Pepticase™	mononitrilase production	395
<i>Rhodotorulla pilimanae</i>	Hy-Soy™	growth medium	628
<i>Saccharomonospora spp.</i>	N-Z-Amine™ A	growth medium	428
<i>Saccharomyces cerevisiae</i>	Hy-Case™ Amino	culture medium	604
	Hy-Soy™	growth medium	644, 650, 654
	Hy-Soy™	Production of papillomavirus vaccines	627
	Hy-Yest™ 412	growth medium for plasmid host cells	653
	N-Z-Amine™	genetic studies	391
	N-Z-Amine™	steroid production	19
	N-Z-Case™	Production of fusion protein MGF and IL-3	573
<i>Saccharopolyspora sp.</i>	N-Z-Case™	Insecticide and miticide production	599
<i>Saccharotrix sp.</i>	Amicase™	growth medium, antibiotic production	620
	Hy-Soy™	growth medium, antibiotic production	620
	N-Z-Amine™ A	growth medium, antibiotic production	620
<i>Salmonella choleraesuis</i>	Edamin™ S	stabilization during lyophilization	656
	N-Z-Amine™ AS	stabilization during lyophilization	656
<i>Salmonella newington</i>	N-Z-Amine™	growth medium	81
<i>Salmonella sp.</i>	Pepticase™	diagnostic medium	551
<i>Salmonella typhi</i>	N-Z-Amine™ A	growth medium	81
<i>Salmonella typhimurium</i>	N-Z-Amine™ A	growth medium	81
<i>Sebekia sp.</i>	N-Z-Amine™ A	Cyclosporin production	612
<i>Selenomonas spp.</i>	Pepticase™	growth medium	423
<i>Serratia marcescens</i>	N-Z-Case™	growth medium	84, 87, 104
<i>Shigella boydii</i>	N-Z-Amine™ A	production of colicins	230

List of microorganisms *continued*

Organism	Product	Application	References
<i>Shigella</i> <i>dispar</i>	N-Z-Amine™ A	production of colicins	230
<i>Stachybotrys</i> sp.	N-Z-Amine™	production of HIV protease inhibitors	586
<i>Staphylococcus aureus</i>	Edamin™ S N-Z-Amine™ N-Z-Amine™ A N-Z-Amine™ A	toxin production synthesis of enterotoxins enterotoxin H production enterotoxin production	105 88, 113, 198, 221 600 58, 61-62, 94, 106, 120, 137, 286, 315, 440, 470, 667 439, 510-511, 647
<i>Staphylococcus</i> spp.	N-Z-Amine™ A Edamin™ S N-Z-Amine™ A	growth medium enterotoxin production enterotoxin production	54 65, 74, 281
<i>Streptococcus agalactiae</i>	Hy-Soy™	growth medium	304
<i>Streptococcus cremoris</i>	Edamin™ S N-Z-Amine™	lactic acid production starter media for cheese cultures	103 468
<i>Streptococcus diacetilactis</i>	N-Z-Amine™ A	lactic acid production	103
<i>Streptococcus equisimilis</i>	N-Z-Case™	synthesis of glutamate dehydrogenase	348
<i>Streptococcus faecalis</i>	N-Z-Amine™ A	streptokinase production	557
<i>Streptococcus mutans</i>	N-Z-Case™	enzyme production	57, 122, 348
<i>Streptococcus sobrinus</i>	N-Z-Case™	growth medium	335
<i>Streptococcus</i> spp	Hy-Case™ SF Hy-Case™ SF N-Z-Amine™ N-Z-Amine™ A	antigen production for vaccines growth medium production of streptokinase, streptodornase antibiotic production	339 204 29 623
<i>Streptomyces albidoflavus</i>	N-Z-Amine™ A	production of biofertilizer	438
<i>Streptomyces alboniger</i>	Hy-Case™ SF	growth medium	258
<i>Streptomyces albus</i>	N-Z-Amine™ A	antibiotic production, albonoursin production	191, 302, 341, 636
<i>Streptomyces ambofaciens</i>	N-Z-Amine™ A	production of spiramycin	472
<i>Streptomyces antibioticus</i>	Hy-Soy™ T N-Z-Amine™ N-Z-Amine™ A N-Z-Amine™ A	antibiotic production synthesis of actinomycin production of antibiotics growth medium	185 98, 132, 427 185, 485 42, 75, 211, 437, 519
<i>Streptomyces aureofaciens</i>	Hy-Case™ Amino Hy-Soy™ N-Z-Amine™ A	production of narasin production of narasin antibiotic production	153 153 141, 153, 164, 168, 179-180, 192
<i>Streptomyces avermitilis</i>	Hy-Case™ SF N-Z-Amine™ N-Z-Amine™ A Peptonized Milk	growth medium production of ivermectin synthesis of ivermectin growth medium	466 388 322 466
<i>Streptomyces avidinii</i>	N-Z-Amine™ A	production of streptavidin	501
<i>Streptomyces azureus</i>	N-Z-Amine™ A	growth medium	220, 351
<i>Streptomyces cattleya</i>	N-Z-Amine™ A	synthesis of thienamycin	253
<i>Streptomyces chromofuscus</i>	N-Z-Amine™ A	production of angiotensin	280
<i>Streptomyces clavuligerus</i>	N-Z-Amine™ A N-Z-Amine™ A	cephalosporin production clavulanic acid production	90, 353 619
<i>Streptomyces coelicolor</i>	Hy-Case™ SF	inoculum medium	665
<i>Streptomyces endus</i>	N-Z-Amine™ A	antibiotic production	252, 387
<i>Streptomyces erythreus</i>	N-Z-Amine™ A	synthesis of erythromycin	23
<i>Streptomyces fradiae</i>	N-Z-Amine™ A N-Z-Amine™ A	synthesis of neomycin growth medium	14, 17 437
<i>Streptomyces gardneri</i>	Hy-Case™ Amino	production of antibiotics	543
<i>Streptomyces graminofaciens</i>	Hy-Soy™	production of streptogramin	16
<i>Streptomyces griseorubiginosus</i>	N-Z-Case™	production of antibiotics	486
<i>Streptomyces griseus</i>	Hy-Case™ SF N-Z-Amine™ N-Z-Amine™ N-Z-Amine™ A	growth medium sporulation medium antibiotic production growth medium	258 465 91, 169 150, 306
<i>Streptomyces halstedii</i>	N-Z-Amine™ A	antibiotic production	278, 338
<i>Streptomyces humidus</i>	N-Z-Amine™ A	antibiotic production	243
<i>Streptomyces hygroscopicus</i>	Hy-Case™ SF	fermentation medium	539
<i>Streptomyces karnatakensis</i>	N-Z-Amine™ A	growth medium	598
<i>Streptomyces kasugaensis</i>	N-Z-Amine™ A	synthesis of kasugamycin	247

<i>Streptomyces lavendulae</i>	N-Z-Amine™	synthesis of moranoline for blood glucose control	349
<i>Streptomyces lipmanii</i>	N-Z-Amine™ A	antibiotic production	145
<i>Streptomyces lividans</i>	N-Z-Amine™ A	genetic studies	292
<i>Streptomyces lomondensis</i>	N-Z-Case™	antibiotic production	400
<i>Streptomyces longwoodensis</i>	N-Z-Amine™ A	growth medium	392
<i>Streptomyces luridus</i>	N-Z-Amine™ A	antibiotic production	342
<i>Streptomyces malachitofucus</i>	Hy-Soy™ T	antibiotic production	224
<i>Streptomyces olivaceus</i>	Primatone™ HS	growth medium	515
<i>Streptomyces parvulus</i>	N-Z-Amine™	growth medium	259, 444, 474
<i>Streptomyces parvus</i>	Hy-Soy™	production of fungicide, eulicin	48
<i>Streptomyces plicatus</i>	N-Z-Amine™ A	antibiotic production	663
<i>Streptomyces pulveraceus</i>	N-Z-Amine™ A	antibiotic production	381
	N-Z-Case™	antibiotic production	381
<i>Streptomyces rochei</i>	N-Z-Amine™	enzyme production	607
<i>Streptomyces roseochromogenus</i>	N-Z-Amine™ A	production of albocycline	248
<i>Streptomyces roseosporus</i>	N-Z-Amine™ A	antibiotic production	469
<i>Streptomyces salvialis</i>	N-Z-Amine™ A	antibiotic production	566
<i>Streptomyces sp</i>	Hy-Case™ SF	production of angiotensin-II antagonist	537, 549, 560-561
	Hy-Case™ SF	fermentation medium	559
	Hy-Soy™	synthesis of antibiotics	365
	N-Z-Amine™	production of angiotensin, converting enzyme inhibitor	319
	N-Z-Amine™	production of ACE inhibitor	318
	N-Z-Amine™	production of elastase inhibitors	571
	N-Z-Amine™	production of xylanase	347
	N-Z-Amine™ A	production of glucuronides	623
	N-Z-Amine™ A	production of nisamycin	564
	N-Z-Amine™ A	antibiotic production	46, 195, 318, 320, 325, 368, 370, 389, 424, 555, 568, 574
			513
<i>Streptomyces tenebrarius</i>	Primatone™ HS	antibiotic production	139, 201
<i>Streptomyces tenjamariensis</i>	N-Z-Amine™ A	production of nebramycin	207
<i>Streptomyces threomyceticus</i>	N-Z-Amine™	growth medium	249
<i>Streptomyces toyocaensis</i>	N-Z-Amine™ A	threomycin production	386
<i>Streptomyces virginiae</i>	N-Z-Amine™ A	antibiotic production	364, 411
	N-Z-Amine™ A	growth medium	435
<i>Streptomyces viridochromogenes</i>	Hy-Case™	growth medium	142
	N-Z-Amine™ A	sulfomycin production	250
<i>Streptomyces viridodiasstaticus</i>	N-Z-Amine™ A	antibiotic production	579, 590
<i>Streptosporangium fragile</i>	N-Z-Amine™ A	synthesis of fragilomycin	299
<i>Streptosporangium spp.</i>	N-Z-Amine™ A	growth medium	428
<i>Streptoverticillium spp.</i>	N-Z-Amine™ A	growth medium	428
<i>Syncephalus nodosa</i>	Edamin™	oxygenation of steroids	15
<i>Tetrahymena pyriformis</i>	Hy-Case™ SF	assay for riboflavin	60, 70
<i>Thermoactinomyces spp.</i>	N-Z-Amine™ A	growth medium	69, 428
<i>Thermomonospora spp.</i>	N-Z-Amine™ A	growth medium	428
<i>Tilletta caries</i>	Edamin™	growth medium	116
<i>Tilletta controversa</i>	Edamin™	growth medium	116
<i>Trichoderma reesei</i>	N-Z-Amine™ A	growth medium	392
<i>Trichoderma viride</i>	Edamin™	hydroxylation of steroids	32
<i>Trypanosoma mega</i>	N-Z-Amine™ AS	growth medium	117
<i>Veillonella alcalescens</i>	N-Z-Amine™ A	growth medium	398
<i>Verticillium sp.</i>	Peptidase™	growth medium	670
<i>Vibrio cholerae</i>	Hy-Case™ SF	growth medium	255, 520
<i>Vibrio harveyi</i>	N-Z-Amine™	genetic studies	459
<i>Xanthomonas campestris</i>	Hy-Soy™ T	production of xanthan gum	181, 238, 242
	N-Z-Amine™ A	production of xanthan gum	297
	N-Z-Amine™ YT	production of xanthan gum	181, 238, 242
<i>Yersinia pestis</i>	N-Z-Amine™	growth medium	209, 263, 330
	N-Z-Amine™ A	growth medium	287, 317

Cell line references

CHO

- Babcock J, Smith S, Huttinga H & Merrill DA. (2007). "Enhancing performance in cell culture: Hydrolysates from plants for medium supplementation." *Genetic Engineering & Biotechnology News* **27**(20): 47-48. (**CHO-K1-origin**)
- Babcock JF, Merrill DA, Smith SR. (2010). A novel approach to the production of plant-derived hydrolysates yields medium supplements with enhanced performance in cell culture systems. *Cells and Culture*. T. Noll, Springer Netherlands. Chapter IV, pp. 847-852. (**CHO-K1-origin**)
- Ballez JS, Mols J, Burteau C, Agathos SN, Schneider YJ. (2004). "Plant protein hydrolysates support CHO-320 cell proliferation and recombinant IFN- γ production in suspension and inside microcarriers in protein-free media." *Cytotechnology* **44**(3):103. (**CHO-K1-origin**)
- Burteau CC, Verhoeve F, Mols JF, Ballez JS, Agathos SN, Schneider YJ. (2003). "Fortification of a protein-free cell culture medium with plant peptones improves cultivation and productivity of an interferon-gamma-producing CHO cell line." *In Vitro Cell and Developmental Biology. Animal.* **39**(7): 291-296. (**CHO-K1-origin**)
- Chottea V, Wåhlberg C, Pettersson, H. (2007). Effect of peptones and study of feeding strategies in a CHO based fed-batch process for the production of a human monoclonal antibody. *Cell Technology for Cell Products*. R. Smith, Springer Netherlands: 371-374. (**CHO-unspecified**)
- Gu X, Xie L, Harmon BJ, Wang DIC. (1997). "Influence of Primate RL supplementation on sialylation of recombinant human interferon- γ produced by Chinese hamster ovary cell culture using serum-free media." *Biotechnology and Bioengineering* **56**(4): 353-360. (**CHO-dhfr--unspecified**)
- Kim DY, Lee J C, Chang HN, Oha DJ (2006). "Development of serum-free media for a recombinant CHO cell line producing recombinant antibody." *Enzyme and Microbial Technology* **39**(3): 426-433. (**CHO-DG44-origin**)
- Kim S. and Lee G. (2009). "Development of serum-free medium supplemented with hydrolysates for the production of therapeutic antibodies in CHO cell cultures using design of experiments." *Applied Microbiology and Biotechnology* **83**(4): 639-648. (**CHO-DG44-origin**)
- Michiels J F, Barbau J, De Boel S, Dessim S, Agathos SN, Schneider YJ. (2011). "Characterisation of beneficial and detrimental effects of a soy peptone, as an additive for CHO cell cultivation." *Process Biochemistry* **46**(3): 671-681. (**CHO-320-origin**)
- Sung YH, Lim SW, Chung JY, Lee GM (2004). "Yeast hydrolysate as a low-cost additive to serum-free medium for the production of human thrombopoietin in suspension cultures of Chinese hamster ovary cells." *Applied Microbiology and Biotechnology* **63**(5): 527-536. (**CHO-DUX-B11-origin**)

Hybridoma

- Babcock JF, Antosh A. (2010a). "Performance enhancing synergy between a wheat hydrolysate and recombinant human serum albumin in SP2/0 hybridoma Cells." *BioProcess International* **8**(7):18. (**SP2/0**)
- Bonarius HPJ, Hatzimanikatis V, Meesters KPH, de Gooijer CD, Schmid G, and Tramper J. (1996). "Metabolic flux analysis of hybridoma cells in different culture media using mass balances." *Biotechnology and Bioengineering* **50**: 299-318. (**Murine**)
- Ganglberger P, Obermüller B, Kainer M, Hinterleitner P, Doblhoff O, and Landauer K. (2007). Optimization of culture medium with the use of protein hydrolysates. In *Cell Technology for Cell Products*. P. Smith, ed. Springer Netherlands. Chapter VI, pp. 553-557. (**SP2/0**)
- Iding, K., Büntemeyer H, Gundermann, F, Deutschmann SM, Kionka C, and Lehmann J. (2001). "An automatic system for the assessment of complex medium additives under cultivation conditions." *Biotechnology and Bioengineering* **73**(6): 442. (**Murine**)
- Jan D, Jones SJ, Emery AN, Al-Rubeai M. (1994). "Peptone, a low-cost growth-promoting nutrient for intensive animal cell culture." *Cytotechnology* **16**(1): 17. (**NS1-derived**)

Schlaeger EJ. (1996). "The protein hydrolysate, Primatone RL, is a cost-effective multiple growth promoter of mammalian cell culture in serum-containing and serum-free media and displays anti-apoptosis properties." *Journal of Immunological Methods* **194**(2): 191. (**Murine-various**)

Zhang Y, Zhou Y, Yu J. (1994). "Effects of peptone on hybridoma growth and monoclonal antibody formation." *Cytotechnology* **16**(3): 147-150. (**Murine-WuT3**)

Vaccine lines

Danlybaeva G, Mazurkova N, Matyushina R, Troshkova G and Podchernyaeva R. (2009). "Nutrient medium from rice flour hydrolysate for cell culturing and evaluation of the sensitivity to influenza viruses." *Bulletin of Experimental Biology and Medicine* **147**(4): 547-550. (**MDCK, Vero**)

Heidemann R, Zhang C, Qi H, Rule J L, et al. (2000). "The use of peptones as medium additives for the production of a recombinant therapeutic protein in high density perfusion cultures of mammalian cells." *Cytotechnology* **32**(2): 157. (**BHK**)

Keay L. (1975). "Autoclavable low cost serum free cell culture media: The growth of L cells and BHK cells on peptones." *Biotechnology and Bioengineering* **17**(5): 745. (**L-cells, BHK**)

Keay L. (1977). "The growth of L-cells and Vero cells on an autoclavable MEM-peptone medium." *Biotechnology and Bioengineering* **19**(3): 399-411. (**L-cells, Vero**)

Mazurkova NA, Kolokol'tsova TD, Nechaeva EA, Shishkina LN, Sergeev AN. (2008). "The Use of components of plant origin in the development of production technology for live cold-adapted cultural influenza vaccine." *Bulletin of Experimental Biology Medicine* **146**(1): 144-147. (**MDCK**)

Mazurkova N, Troshkova G, Sumkina, T, Kolokol'tsova, T, et al. (2008a). "Comparative analysis of reproduction of influenza virus strains in cell lines perspective for the creation of cultural vaccines grown on nutrient medium on the basis of rice flour protein hydrolysate." *Bulletin of Experimental Biology and Medicine* **146**(4): 547-550. (**MDCK, Vero**)

Mazurkova N, Ryabchikova E, Troshkova G, Getmanova T, et al. (2009). "Morphological and proliferative characteristics of Vero and MDCK cells during culturing in nutrient media on the basis of hydrolysates of plant proteins." *Bulletin of Experimental Biology and Medicine* **147**(4): 551-554. (**MDCK, Vero**)

Rourou, S., A. van der Ark, et al. (2009). "Development of an animal-component free medium for Vero cell culture." *Biotechnology Progress* **25**(6): 1752. (**Vero**)

Saha, S. N. and A. K. Sen (1989). "Studies on the development of a medium with peptone and casein hydrolysate for the production of foot-and-mouth disease vaccine in BHK-21 cells." *Vaccine* **7**(4): 357. (**BHK**)

Saha, S. N. and A. K. Sen (1989a). "Partial replacement of serum with peptone and lactalbumin hydrolysate for the production of foot-and-mouth disease vaccine in BHK-21 cells." *Acta Virologica* **33**(4): 338. (**BHK**)

Human

Chun BH, Lee YK, Bang WG, Chung N. (2005). "Use of plant protein hydrolysates for varicella virus production in serum-free medium." *Biotechnology Letters* **27**(4): 243. (**MRC-5**)

Han X, Sun L, Fang Q, Li D, Gong X, Wu Y, Yang S, and Bing S. (2007). "Transient expression of osteopontin in HEK 293 cells in serum-free culture." *Enzyme and microbial technology* **41**: 133-140. (**HEK 293**)

Lee J, Lee J, Hwang H, Jung E, Huh S, Hyun J and Park D. (2009). "Promotion of stem cell proliferation by vegetable peptone." *Cell Proliferation* **42**(5): 595. (**Human CB-MSC, Human ADSC**)

Pham PL, Perret S, Cass B, Carpentier E, St.-Laurent, G, et al. (2003). "Large-scale transient transfection of serum-free suspension-growing HEK293 EBNA1 cells: Peptone additives improve cell growth and transfection efficiency." *Biotechnology and Bioengineering* **84**(3): 332-342. (**HEK 293**)

Cell line references *continued*

Pham PL, Perret S, Doan HC, Cass B, St-Laurent G, et al. (2005). "Transient gene expression in HEK293 cells: Peptone addition post-transfection improves recombinant protein synthesis." *Biotechnology and Bioengineering* **90**(3): 332. (HEK 293)

Insect cells

Agathos SN. (2007). Development of serum-free media for lepidopteran insect cell lines. *Methods in Molecular Biology*. **388**: 155. (*Spodoptera frugiperda* and *Trichoplusia ni*)

Donaldson MS, Shuler ML. (1998). "Low-cost serum-free medium for the BTI-Tn5B1-4 insect cell line." *Biotechnology Progress* **14**(4): 573-579. (BTI-Tn5B1-4 [High-Five])

Galesi A, Pereira CA, Moraes AM (2007). "Culture of transgenic Drosophila melanogaster Schneider 2 cells in serum-free media based on TC100 basal medium." *Biotechnology Journal* **2**(11): 1399-1407. (*Drosophila melanogaster*)

Ikonomou L, Bastin G, Schneider YJ, Agathos SN (2001). "Design of an efficient medium for insect cell growth and recombinant protein production." *In Vitro Cell and Developmental Biology. Animal*. **37**: 549-559. (Sf-9 and High-Five)

Ikonomou, L, Schneider YJ, Agathos SN. (2003). "Insect cell culture for industrial production of recombinant proteins." *Applied Microbiology and Biotechnology* **62**(1): 1-20. (Various-review)

Kwon MS, D. T., Park EY (2005). "Use of plant-derived protein hydrolysates for enhancing growth of Bombyx mori insect cells in suspension culture." *Biotechnology and Applied Biochemistry*. (*Bombyx mori*)

Marteijn RCL, Jurrius O, Dhont J, de Gooijer CD, Tramper J, and Martens DE. (2003). "Optimization of a feed medium for fed-batch culture of insect cells using a genetic algorithm." *Biotechnology and Bioengineering* **81**(3): 269-278. (Sf-9)

Mendonça RC, De Oliveira EC, Pereira CA, Lebrun I (2007). "Effect of bioactive peptides isolated from yeastolate, lactalbumin and N-Z-Case in the insect cell growth." *Bioprocess and Biosystems Engineering* **30**(3): 157-164. (Sf-9)

Park EY, Kageshima A, Kwon MS, Kato T. (2007). "Enhanced production of secretory β 1,3-N-acetylglucosaminyltransferase 2 fusion protein into hemolymph of *Bombyx mori* larvae using recombinant BmNPV bacmid integrated signal sequence." *Journal of Biotechnology* **129**(4): 681-688. (*Bombyx mori*)

Plant cells

Anjum S, Zia M, Chaudhary MF. (2006). "Investigations of different strategies for high frequency regeneration of *Dendrobium malones* 'Victory.'" *African Journal of Biotechnology* **5**(19): 1738-1743. (*Dendrobium malones*)

Gamborg OL, Larue TAG. (1971). "Ethylene Production by Plant Cell Cultures: The Effect of Auxins, Abscisic Acid, and Kinetin on Ethylene Production in Suspension Cultures of Rose and Ruta Cells." *Plant Physiology* **48**(4): 399-401. (*Rosa sp.* and *Ruta sp.*)

Gamborg OL, Miller RA, Ojima K. (1968). "Nutrient requirements of suspension cultures of soybean root cells." *Experimental Cell Research* **50**(1): 151. (*Glycine max*)

Parc G, Rembur J, Rech P, Chiriqui D. (2007). "In vitro culture of tobacco callus on medium containing peptone and phytate leads to growth improvement and higher genetic stability." *Plant Cell Reports* **26**(2): 145. (*Nicotiana tabacum*)

Ranch JP, Giles KL. (1980). "Factors Affecting Growth and Aggregate Dissociation in Batch Suspension Cultures of *Datura innoxia* (Miller)." *Annals of Botany* **46**(6): 667. (*Datura innoxia*)

References

- 1 Wyman L. and Bowen H.E., 55th Ann. Mtg., Soc. Am. Bacteriol., 1900.
- 2 Mazurs C. et al, ABSTRACT, 1901.
- 3 Vosper N., Appl. Reports in Mammalian Cell Culture, New Brunswick Scientific Co., Inc. CC106. 1902.
- 4 Vosper N., Appl. Reports in Mammalian Cell Culture, New Brunswick Scientific Co., Inc., CC105. 1903.
- 5 Dubos R.J. and Davis B.D., J. Exptl. Med., 1946, 5, 409.
- 6 Dubos R.J. and Middlebrook G., Am. Rev. Tuberculosis, 1947, 56, 334.
- 7 Mueller J.H and Miller P. A., J. Bacteriol., 1948, 56, 219.
- 8 Spies J.R. and Chambers D.C. Anal. Chem., 1949, 21, 1249.
- 9 Eggerth A. H., Proc. Soc. Exptl. Biol. Med. 1950, 73, 542.
- 10 Sabin B.A. et al, US patent 2,516,080, 1950.
- 11 GB 0.650.087, 1951.
- 12 Levine L. and Stone J.L., J. Immunol., 1951, 67, 235.
- 13 GB 0.677.423, 1952.
- 14 Hickey R.J. and Tresner H.D., J. Bacteriol., 1952, 64, 891.
- 15 Murray H.C. et al, US patent 2,602,769. 1952.
- 16 Charney J. et al, Antibiot. and Chemother, 1953, 3, 1283.
- 17 Dulmage H.T., Appl. Microbiol., 1953, 1, 103.
- 18 Mueller J.H and Miller P. A., J. Bacteriol. 1953, 67, 271.
- 19 Dulaney E.L. et al, Appl. Microbiol., 1954, 2, 371.
- 20 Sher I.H. and Mallette M.F., Arch. of Biochem. and Biophys., 1954, 52, 331.
- 21 Sher I.H. and Mallette M.F., Arch. of Biochem. and Biophys., 1954, 53, 354.
- 22 Sher I.H. and Mallette M.F., Arch. of Biochem. and Biophys., 1954, 53, 370.
- 23 Pettinga C.W. et al, J. Am. Chem. Soc., 1954, 76, 569.
- 24 Tanner F.W. et al, US patent 2,691,618., 1954.
- 25 Dulaney E.L. et al, Appl. Microbiol., 1955, 3, 336.
- 26 Gyorgy P. and Rose C.S., J. Bacteriol., 1955, 69, 483.
- 27 Mueller J.H and Miller P. A., J. Bacteriol., 1955, 69, 634.
- 28 Bachrach H.L. et al, Science, 1955, 122, 1269.
- 29 Ablondi F.B. and Adam J.N. Jr., US patent 2,701,227, 1955.
- 30 Rickes E.L. and Wood T.R., US patent 2,703,302, 1955.
- 31 Shull G.M. et al, US patent 2,721,163, 1955.
- 32 McAleer W.J. and Dulaney E.L., Arch. Biochem. and Biophys., 1956, 62, 111.
- 33 Hotta S. and Evans C.A., J. Infect. Diseases, 1956, 98, 88.
- 34 Lein J. and Gourevitch A., US patent 2,739,924, 1956.
- 35 Sabin B.A. et al, US patent 2,757,123, 1956.
- 36 Orlean et all., Journal of Systematic Bacteriol., 1978, vol 28 nr2/ p194-196.
- 37 Duff J.T. et al, J. Bacteriol., 1957, 42, 73.
- 38 Howe C. et al, J. Bacteriol., 1957, 74, 365.
- 39 Bauman N. and Davis B.D., Science., 1957, 126,170.
- 40 Shantz E.M. and Steward F.C., Ann. Bot., N. S., 1958, 23, 371.
- 41 Cardella M.A. et al, J. Bacteriol., 1958, 75, 360.
- 42 Katz E. and Goss W.A., Biochem. J., 1959, 73, 458.
- 43 Kheshgi S. and Saunders R.J., Biochem. and Microbiol. Technol. and Eng., 1959, 1, 115.
- 44 Shull G.M., US patent 2,890,153, 1959.
- 45 Baldacci D. et al, J. Bacteriol., 1960, 79, 379.
- 46 Gottlieb D., Appl. Microbiol., 1961, 9, 55.
- 47 Lechevalier H.A. et al, J. Gen. Microbiol., 1961, 26, 11.
- 48 Charney J. et al, US patent 2,998,438, 1961.
- 49 Latham W.C. et al, Appl. Microbiol., 1962, 10, 146.
- 50 Maretzki A. and Mallette M.F., J. Bacteriol., 1962, 83, 720.
- 51 Murashige T and Skoog F., Phys. Plant., 1962, 15, 473.
- 52 Luedemann G.M. and Brodsky B.C., Antimicrob. Ag. and Chemother., 1963, 116.
- 53 Pollock M.E. et al, Exp. Cell Res., 1963, 31, 321.
- 54 Casman E.P. and Bennett R. W., J. Bacteriol., 1963, 86, 18.
- 55 Higuchi K., J. Infect. Diseases., 1963, 112, 213.
- 56 Luedemann G.M. and Brodsky B.C., Antimicrob. Ag. and Chemother., 1964, 47.
- 57 Bleiweis A.S. and Zimmerman L.N., J. Bacteriol., 1964, 88, 654.
- 58 Schantz E.J. et al, Biochem., 1965, 4, 1011.
- 59 Bergdoll M.S. et al, J. of Bacteriol., 1965, 90, 1481.
- 60 Baker H. et al, Am. J. Clin Nutr., 1966, 19, 17.
- 61 Kato E. et al, Appl. Microbiol., 1966, 966.
- 62 Rosenwald A. J. and Lincoln R.E., J. Bacteriol., 1966, 92, 279.
- 63 Huang H.T. and English A. R., US patent 3,239,427, 1966.
- 64 Buskirk H.H., Appl. Microbiol., 1967, 15, 1442.
- 65 Casman E.P. et al, J. Bacteriol., 1967, 94, 1875.
- 66 Sheppard D.E. and Englesberg E., J. Mol. Biol., 1967, 25, 443.
- 67 Stainer D.W. et al, Can. J. Microbiol., 1968, 14, 1155.
- 68 Luedemann G.M., J. Bacteriol., 1968, 96, 1848.
- 69 Lechevalier M.P., J. lab. and Clin. Mad., 1968, 71, 934.
- 70 Cox D. et al, J. Protozool., 1968, 15, 713.
- 71 Zygmunt W.A. and Browder H.P., US patent 3,398,056, 1968.
- 72 Moravek L. et al, J. Biol. Chem., 1969, 244, 497.
- 73 Weinstein M.J. et al, US patent 3,454,696, 1969.
- 74 Bergdoll M.S. et al, Academic Press, New York., 1970, 3, 265.
- 75 Jones G.H. and Weissbach H., Arch. of Biochem. and Biophys., 1970, 137, 558.
- 76 Luedemann G.M. and Weinstein M. J., US patent 3,499,078, 1970.
- 77 Luedemann G.M. and Weinstein M. J., US patent 3,511,909, 1970.
- 78 Fung D.Y.C. and Wagner J., Appl. Microbiol., 1971, 21, 559.
- 79 Clifford W.J. and Anellis A., Appl. Microbiol., 1971, 22, 856.
- 80 Kele R.A. and McCoy E., Appl. Microbiol., 1971, 22, 728.
- 81 Jay J.M, Infect. and Immun., 1971, 3, 544.
- 82 Gruber J., Infect. and Immun., 1971, 3, 574.
- 83 Luedemann G.M., Int. J. System. Bacteriol., 1971, 21, 240.
- 84 Williams R.P. et al, J. Bacteriol., 1971, 106, 444.
- 85 Das Gupta B.A., J. Bacteriol., 1971, 108, 1051.
- 86 Coats J.H. and Argoudelis A.D., J. Bacteriol., 1971, 108, 459.
- 87 Williams R.P. et al, J. Bacteriol., 1971, 438.
- 88 Hurst A. and Kruse H., Ag. and Chemother., 1972, 1, 277.
- 89 Kominek L.A., Antimicrob. Ag. and Chemother., 1972, 1, 123.
- 90 J. G. Whitney J.G. et al, Antimicrob. Ag. and Chemother., 1972, 1, 247.
- 91 Stapley E.O. et al, Antimicrob. Ag. and Chemother., 1972, 2, 122.
- 92 Peterson R.E. and Ciegler A., Appl. Microbiol., 1972, 23, 671.
- 93 Marshek W.J. et al, Appl. Microbiol., 1972, 23, 72.
- 94 Dietrich G.G. et al, Appl. Microbiol., 1972, 24, 561.
- 95 Schantz E.J. et al, Biochem., 1972, 11, 360.
- 96 Tomlinson G.A. and Hochstein L.I., Can. J. Microbiol., 1972, 18, 1973,
- 97 G. H. Wagman G.H. et al, J. Antibiot., 1972, 25, 641.
- 98 Gallo M. and Katz E., J. Bacteriol., 1972, 109, 659.
- 99 Eleuterio M. et al, J. Bacteriol., 1972, 111, 383.
- 100 Klein S., J. Protozool., 1972, 19, 140.
- 101 Joshi P.C. and Johri B.M., Phytomorphology., 1972, 22, 195.
- 102 Weinstein M.J. et al, US patent 3,632,750., 1972.
- 103 Molskness T.A. et al, Appl. Microbiol., 1973, 25, 373.
- 104 Williams R.P., Appl. Microbiol., 1973, 25, 396.
- 105 Vandenbosch L.L. et al, Appl. Microbiol., 1973, 25, 498.
- 106 Joseph F. Metzger J.F. et al, Appl. Microbiol., 1973, 25, 770.
- 107 Miller R.D and Fung D.Y.C., Appl. Microbiol., 1973, 25, 800.
- 108 Fung D.Y.C., Appl. Microbiol., 1973, 26, 638.
- 109 Cascieri T. Jr. and Mallette M.F., Appl. Microbiol., 1973, 26, 975.
- 110 Woodburn M. et al, Appl. Microbiol., 1973, 25, 825.
- 111 Stainer D.W. and Scholte J.M.C., Biotech. & Bioeng. Symp. No. 4., 1973, 283.
- 112 Lasfargues E.Y. et al, In Vitro, 1973, 8, 494.
- 113 Hurst A. et al, J. Bacteriol., 1973, 116, 901.
- 114 Butenko R.G. et al, Phytomorphology., 1973, 22, 140.

References *continued*

- 115 Baker H. et al, Am. J. Clin. Nutr., 1974, 21, 322.
116 Trione E.J., Amer. J. Bot., 1974, 61, 914.
117 Bacchi C.J. et al, Antimicrob. Ag. and Chemother., 1974, 6, 785.
118 Wagman G.H. et al, Antimicrob. Ag. and Chemother., 1974, 6, 144.
119 Barksdale A.W. and Lasure L.L., Appl. Microbiol., 1974, 28, 544.
120 Carpenter D.F. and Silverman G.J., Appl. Microbiol., 1974, 28, 628.
121 Robbins R. et al, Appl. Microbiol., 1974, 28, 946.
122 Appelbaum B. and Zimmerman L.N., Infect. and Immun., 1974, 10, 991.
123 Barksdale A.W. et al, J. Gen. Microbiol., 1974, 82, 295.
124 Goldberg B. et al, J. Protozool., 1974, 21, 322.
125 Calnek B.W. and Hitchner S. B., US patent 3,783,098., 1974.
126 Merkel J.R. et al, Appl. Microbiol., 1975, 29, 145.
127 Yang K.H. and Sugiyama H., Appl. Microbiol., 1975, 29, 598.
128 Keay L., Biotech. and Bioeng., 1975, 17, 745.
129 Nelson G.E.N. et al, Dev. in Ind. Microbiol., 1975, 16, 385.
130 Andron L.A. and Eigelsbach H.T., Infect. and Immun., 1975, 12, 137.
131 Metzger J.F. et al, Infect. and Immun., 1975, 12, 93.
132 Jones G. H., J. Bacteriol., 1975, 124, 364.
133 Coats J.H., Methods in Enzymology., 1975, 43, 755.
134 Celmer W.D. et al, US patent 3,884,763, 1975.
135 Celmer W.D. et al, US patent 3,914,218, 1975.
136 Jones G. H., Biochem., 1976, 15, 3331.
137 Carpenter D.F. et al, Dev. in Ind. Microbiol., 1976, 17, 363.
138 Nelson G.E.N. and Peterson R.E., Dev. in Ind. Microbiol., 1976, 17, 399.
139 Stark W.M. et al, Dev. in Ind. Microbiol., 1976, 17, 61.
140 GB 1,454,735, 1976.
141 Dorman D.E. et al, Heli. Chem., Acta., 1976, 59, 2625.
142 Hirsch C.F. and Ensign J.C., J. Bacteriol., 1976, 126, 13.
143 Nara T. et al, US patent 3,931,400, 1976.
144 Daum S.J. and Clarke R.L., US patent 3,972,930, 1976.
145 Hamill R.L. et al, US patent 3,973,015, 1976.
146 Daum S.J. and Clarke R.L., US patent 3,982,996, 1976.
147 Berg D.H. et al, US patent 3,991,052, 1976.
148 Celmer W.D. et al, US patent 3,991,183, 1976.
149 Ort D.S., Appl. and Environ. Microbiol., 1977, 34, 710.
150 Gerber N.N. and Lechevalier H.A., Appl. and Environ. Microbiol., 1977, 34, 857.
151 Mizrahi A., Biotech. and Bioeng., 1977, 19, 1557.
152 Keay L., Biotech. and Bioeng., 1977, 19, 399.
153 Boeck L.D. et al, Dev. in Ind. Microbiol., 1977, 18, 471.
154 Nelson G.E.N. and Peterson R.E., Dev. in Ind. Microbiol., 1977, 8, 511.
155 Hink W.F. et al, In Vitro., 1977, 13, 177.
156 Landolo J.J. and Shafer W.M., Infect. and Immun., 1977, 16, 610.
157 Keay L. and Hertelendy F., IRCS Med. Sci.: Cancer; Cell and Membrane Biol.; Endocrine Sys., 1977, 5, 176.
158 Argoudelis A. D. et al, J. Antibiot., 1977, 30, 474.
159 Mizrahi A. and Avihoo A., J. Biol. Std., 1977, 5, 31.
160 N. Ausavanodom N. et al, J. Dairy Sci., 1977, 60, 1245.
161 Blattner F.R. et al, Science, 1977, 196, 161.
162 Berg D.H. et al, US patent 4,001,086, 1977.
163 Chaiet L. et al, US patent 4,007,090, 1977.
164 Boeck L.D. and Berg D.H., US patent 4,009,262, 1977.
165 Hoehn M.M and Michel K. H., US patent 4,024,245, 1977.
166 Daum S.L. and Clarke R.L., US patent 4,028,188, 1977.
167 Celmer W.P. et al, US patent 4,032,631, 1977.
168 Berg D.H. et al, US patent 4,035,481, 1977.
169 Stapley E.O. and Mata J. M., US patent 4,036,696, 1977.
170 Nara T. et al, US patent 4,045,298, 1977.
171 Shafer W.M. and Landolo J.J., Appl. and Environ. Microbiol., 1978, 36, 389.
172 Merkel J.R. and Dreisbach J. H., Biochem., 1978, 17, 2857.
173 Nelson G.E.N. and Peterson R.E., Dev. in Ind. Microbiol., 1978, 19, 451.
174 Shafer W.M. and Landolo J.J., Infect. and Immun., 1978, 20, 273.
175 Hasegawa T. et al, Int. J. Syst. Bacteriol., 1978, 28, 304.
176 Dreisbach J.H. and Merkel J.R., J. Bacteriol., 1978, 135, 521.
177 Anonymous, Procedures for Use of Charon Phages in Recombinant DNA Molecule Research, National Institute of Allergy and Infectious Diseases, NIH, Bethesda, MD., 1978.
178 Blattner F.R. et al, Science, 1978, 202, 1279.
179 Nakatsukasa W.M. and Hamill R.L., US patent 4,110,435, 1978.
180 Nakatsukasa W.M. and Hamill R.L., US patent 4,110,436, 1978.
181 Werneau W.C., US patent 4,119,546, 1978.
182 Burg R.W. et al, Antimicrob. Ag. and Chemother., 1979, 15, 361.
183 Spira W.M. and Silverman G.J., Appl. and Environ. Microbiol., 1979, 37, 109.
184 Walsh R.G. et. al., Int. J. Pept. Prot. Res., 1979, 14, 290.
185 Lui C. et al, J. Antibiot., 1979, 32, 95.
186 Robinson J.M. et al, J. Bacteriol., 1979, 137, 1158.
187 Zoon K.C. et al, J. Gen. Virol., 1979, 44, 227.
188 Williams B.G. and Blattner F.R., J. Virol., 1979, 29, 555.
189 Hohn B., Methods in Enzymology, 1979, 68, 299.
190 Jenkin H.M. et al, Proc. Soc. Exptl. Biol. and Med., 1979, 160, 59.
191 Hamill R.L. and Hoehn M.M., US patent 4,133,876, 1979.
192 Nakatsukasa W.M. et al, US patent 4,141,907, 1979.
193 Chabala J.C. et al, US patent 4,171,314, 1979.
194 Chabala J.C. et al, US patent 4,173,571, 1979.
195 Godfrey O.W. et al, US patent 4,180,564, 1979.
196 Ranch J.P. and Giles K.L., Ann. Bot., 1980, 46, 667.
197 Gerber N.N. et al, Antimicrob. Ag. and Chemother., 1980, 17, 636.
198 Pattee P.A. and Glatz B.A., Appl. and Environ. Microbiol., 1980, 39, 186.
199 Mizrahi A. et al, Biotechnol. Lett., 1980, 2, 267.
200 Mafthews D.M. and Payne J.W., Current Topics in Membranes and Transport., 1980, 14, 331.
201 Stark W.M. et al, Dev. in Ind. Microbiol., 1980, 21, 77.
202 Reuveny S. et al, Develop. Biol. Std., 1980, 46, 281.
203 Taylor A.F. et al, Gene, 1980, 9, 321.
204 Carey R.B. et al, Infect. and Imm., 1980, 28, 195.
205 Ohkuma H. et al, J. Antibiot., 1980, 33, 1087.
206 Liu W. et al, J. Antibiot., 1980, 33, 1256.
207 Hotta K. et al, J. Antibiot., 1980, 33, 1515.
208 Sato S. et al, J. Antibiot., 1980, 33, 364.
209 Darveau R.P. et al, J. Bacteriol., 1980, 143, 942.
210 Berg D.E. et al, J. Bacteriol., 1980, 142, 439.
211 Jones G.H., J. Biol. Chem., 1980, 255, 1945.
212 Robinson J.M. et al, J. Gen. Microbiol., 1980, 118, 529.
213 Baker H. et al, Methods in Enzymol., 1980, 66, 490.
214 Payne J.W., Microorganisms and Nitrogen Sources, John Wiley and Sons, New York., 1980.
215 Vedamethu E.R., US Patent 4,191,782, 1980.
216 Chabala J.C. and Fisher M.H., US patent 4,199,569, 1980.
217 Fisher M.H. et al, US patent 4,200,581, 980.
218 Mrozik H.H. et al, US patent 4,201,861, 1980.
219 Mann G.F., US patent 4,232,124, 1980.
220 Ogata S. et al, Appl. and Environ. Microbiol., 1981, 42,135.
221 Tweten R.K. and Landolo J.J., Infect. and Imm., 1981, 34, 900.
222 Dyer D.W. and Landolo J.J., Infect. and Immun., 1981, 33, 450.
223 Waitz J.A. et al, J. Antibiot., 1981, 34, 1101.
224 Liu C. et al, J. Antibiot., 1981, 34, 133.
225 Albers-Schonberg G. et al, J. Antibiot., 1981, 34, 507.
226 Autissier D. et al, J. Antibiot., 1981, 34, 536.
227 Mazieres N. et al, J. Antibiot., 1981, 34, 544.
228 Liu W. et al, J. Antibiot., 1981, 34, 791.
229 Comer M.M., J. Bacteriol., 1981, 146, 269.
230 Ferber D.M. et al, J. Bacteriol., 1981, 146, 506.
231 Barnes W.M., J. Bacteriol., 1981, 147, 124.
232 Bergdoll M.S. et al, Lancet., 1981, May, 1017.

- 233 Mizrahi A., Meth. in Enzymol., 1981, 78, 54.
- 234 Sykes R.B. et al, Nature., 1981, 291, 489.
- 235 Phillips G.C. and Collins G.B., Plant Cell Tissue Organ Culture., 1981, 1, 123.
- 236 Sykes R.B. et al, US patent 4,249,008, 1981.
- 237 Pereira DaLuz A., US patent 4,275,155, 1981.
- 238 Werneau W.C., US patent 4,282,321, 1981.
- 239 Yamamoto H., US patent 4,288,547, 1981.
- 240 Boeck L.D. and Kastner R.E., US patent 4,288,549, 1981.
- 241 Albers-Schonberg G. et al, US patent 4,294,846, 1981.
- 242 Werneau W.C., US patent 4,296,203, 1981.
- 243 Hatt H.D., American Type Culture Collection Catalogue of Strains 1, Fifteenth Edition, Rockville, Maryland, 1982, 215.
- 244 Hatt H.D., American Type Culture Collection Catalogue of Strains 1, Fifteenth Edition, Rockville, Maryland, 1982, 129.
- 245 Hatt H.D., American Type Culture Collection Catalogue of Strains 1, Fifteenth Edition, Rockville, Maryland, 1982, 144.
- 246 Hatt H.D., American Type Culture Collection Catalogue of Strains 1, Fifteenth Edition, Rockville, Maryland, 1982, 158.
- 247 Hatt H.D., American Type Culture Collection Catalogue of Strains 1, Fifteenth Edition, Rockville, Maryland, 1982, 216.
- 248 Hatt H.D., American Type Culture Collection Catalogue of Strains 1, Fifteenth Edition, Rockville, Maryland, 1982, 223.
- 249 Hatt H.D., American Type Culture Collection Catalogue of Strains 1, Fifteenth Edition, Rockville, Maryland, 1982, 227.
- 250 Hatt H.D., American Type Culture Collection Catalogue of Strains 1, Fifteenth Edition, Rockville, Maryland, 1982, 229.
- 251 Mallonee D.H. et al, Appl. and Environ. Microbiol., 1982, 43, 397.
- 252 Ogata S. et al, Appl. and Environ. Microbiol., 1982, 43, 1182."
- 253 Foor F. et al, Dev. in Ind. Microbiol., 1982, 23, 305.
- 254 Bloom et al., Appl. and Environ. Microbiol., 1989, Vol. 55, No. 9, p. 2161-2166.
- 255 Cryz Jr. S.J. et al, Infect. and Imm., 1982, 38, 21.
- 256 Smith G.C. and Merkel J.R., Infect. and Immun., 1982, 35, 1155.
- 257 Wells J. S. et al, J. Antibiot., 1982, 35, 184.
- 258 Coleman R.H. and Ensign J.C., J. Bacteriol., 1982, 149, 1102.
- 259 Ochi K., J. Bacteriol., 1982, 150, 598.
- 260 Parkinson J.S. and Houts S.E., J. Bacteriol., 1982, 151, 106.
- 261 Hui I. et al, J. Bacteriol., 1982, 152, 1022.
- 262 Ferrari F.A. et al, J. Bacteriol., 1982, 152, 809.
- 263 Charnetzky W.T. and Brubaker R.E., J. Bacteriol., 1982, 149, 1089.
- 264 Cryz Jr. S.J. et al, J. Clin. Microbiol., 1982, 16, 41.
- 265 Hirschel B.J. and Berg D.E., J. Mol. Biol., 1982, 105, 155.
- 266 Chihiro Sasakawa and Douglas E. Berg., J. Mol. Biol., 1982, 159, 257.
- 267 Motamedi H. et al, J. Mol. Biol., 1982, 162, 535.
- 268 von Meyenburg K. et al, Mol. Gen. Genet., 1982, 188, 240.
- 269 Comer M.M., Mol. Gen. Genet., 1982, 132,187.
- 270 Maniatis T. et al, Molecular Cloning, Cold Spring Harbor Laboratory, Cold Spring Harbor, N.Y., 1982.
- 271 Pickrum H.M., poster presented at 82nd Ann. Mtg. of Am. Soc. Microbiol., 1982.
- 272 Harold D. Hatt, ed. American Type Culture Collection Catalogue of Strains I, Fifteenth, Rockville, MD 1982, 67.
- 273 Albers-Schonberg G. et al, US patent 4,310,519, 1982.
- 274 Hernandez S. et al, US patent 4,311,693, 1982.
- 275 Albers-Schonberg G., US patent 4,319,039, 1982.
- 276 Tejera E. et al, US patent 4,330,624, 1982.
- 277 Albers-Schonberg G. et al, US patent 4,342,767, 1982.
- 278 Celmer W.D. et al, US patent 4,361,649, 1982.
- 279 Gado I. et al, US patent 4,363,876, 1982.
- 280 Wilgus R. et al, Am. Socy. Microbiol., 1983, poster 058.
- 281 Hoover D.G. et al., Appl. and Environ. Microbiol., 1983, 46, 649.
- 282 Gilbert J.P. et al, Appl. and Environ. Microbiol., 1983, 46, 452.
- 283 EP 0.090.886, 1983.
- 284 Berg D.E. et al, Gene, 1983, 1, 22.
- 285 Pang P.P. and Walker G.C., J. Bacteriol., 1983, 153, 1172.
- 286 Tweten R.K. and Landolo J.J., J. Bacteriol., 1983, 153, 297.
- 287 Moncla B.J. et al, J. Bacteriol., 1983, 153, 340.
- 288 Crofton S. and Dennis P.P., J. Bacteriol., 1983, 154, 58.
- 289 Pierce J.R. et al, J. Bacteriol., 1983, 155, 330.
- 290 Slocum M.K. and Parkinson J.S., J. Bacteriol., 1983, 155, 565.
- 291 Golden S.S. and Sherman L.A., J. Bacteriol., 1983, 155, 966.
- 292 Horinouchi S. et al, J. Bacteriol., 1983, 55, 1238.
- 293 Barnes W.M. and Tuley E., J. Mol. Biol., 1983, 165, 443.
- 294 Yamanouchi O.H., JP 08.107.173, 1983.
- 295 Jarai M. et al, US patent 4,369,251, 1983.
- 296 Jarai M. et al, US patent 4,371,622, 1983.
- 297 Richmon J.B., US patent 4,375,512, 1983.
- 298 Albers-Schonberg G.E. et al, US patent 4,385,065, 1983.
- 299 Nash C.H. et al, US patent 4,385,122, 1983.
- 300 Daum S.J. and Clarke R.L., US patent 4,396,609, 1983.
- 301 Marscheck W.J. et al, US patent 4,397,947, 1983.
- 302 Celmer W.D. et al, US patent 4,411,892, 1983.
- 303 Ormond R.E., US patent 4,412,991, 1983.
- 304 Carlo D.J. et al, US patent 4,413,057, 1983.
- 305 Hernandez S., US patent 4,415,669, 1983.
- 306 Nozaki Y. et al, Agric. Biol. Chem., 1984, 48, 37.
- 307 Medwid R.D. and Grant D.W., Appl. and Environ. Microbiol., 1984, 48, 1067.
- 308 Fusee M.C. and Weber J.E., Appl. and Environ. Microbiol., 1984, 48, 694.
- 309 Stein D.C. et al, Appl. and Environ. Microbiol., 1984, 48, 280.
- 310 Anonymous, ATCC Update, Bacteria & Phages, American Type Culture Collection, Rockville, MD., 1984, 4.
- 311 ATCC Update., Bacteria & Phages, American Type Culture Collection, Rockville, MD., 1984, 10.
- 312 Wood S., Bio Techniques., 1984, Nov./Dec., 314.
- 313 Wilson K.E. et al, EP 0.106.408, 1984.
- 314 Kaper J.B. et al, Infect. and Imm., 1984, 45, 290.
- 315 Igarashi H. et al, Infect. and Immun., 1984, 44, 175.
- 316 Reiser R.F. et al, Infect. and Immun., 1984, 45, 625.
- 317 Straley S.C and Hannon P., Infect. and Immun., 1984, 45, 655.
- 318 Lee M.D. et al, J. Antibiot., 1984, 37 (10) 1149-1152.
- 319 Huang L. et al, J. Antibiot., 1984, 37 (5) 462-465.
- 320 Lee M.D. et al, J. Antibiot., 1984, 37, 1149.
- 321 Chaiet L. et al, J. Antibiot., 1984, 37, 207.
- 322 Calcott P.H. and Fatig R.O., J. Antibiot., 1984, 37, 253.
- 323 Patel M. et al, J. Antibiot., 1984, 37, 413.
- 324 Parker W.L. et al, J. Antibiot., 1984, 37, 431.
- 325 Omura S. et al, J. Antibiot., 1984, 37, 700.
- 326 Konishi M. et al, J. Antibiot., 1984, 37, 949.
- 327 Huang L. et al, J. Antibiot., 1984, 37, 970.
- 328 Dykstra C.C. et al, J. Bacteriol., 1984, 157, 21.
- 329 Bruist M.F. and Simon M.L., J. Bacteriol., 1984, 159, 71.
- 330 Goguen J.D. et al, J. Bacteriol., 1984, 60, 842.
- 331 Sutton B.C.S. et al, J. Bacteriol., 1984, 158, 920.
- 332 Cooper J.D.H. et al, J. Chromatogr., 1984, 285, 490.
- 333 Zenk M.H. et al, Naturwissenschaften, 1984, 71 (5) 266.
- 334 Anonymous, Official Methods of Analysis of the Assn. of Off. Anal. Chem. 1984, 14th ed., 2.057.
- 335 Williams S., Official Methods of Analysis of the Association of Official Analytical Chemists., 1984, 42.271.
- 336 Schulte U. et al, Plant Cell Rep., 1984, 3 (2), 51-54.
- 337 Albers-Schonberg G. et al, US patent 4,429,042, 1984.
- 338 Celmer W.D. et al, US patent 4,431,801, 1984.
- 339 Swenson R.M. et al, US patent 4,439,422, 1984.
- 340 Waitz J.A. et al, US patent 4,440,751, 1984.
- 341 Celmer W.D. et al, US patent 4,450,237, 1984.
- 342 Johnson R.D. et al, US patent 4,463,092, 1984.
- 343 Kempf A.J. et al, US patent 4,476,139, 1984.
- 344 Gonzalez C.F., US patent 4,477,471, 1984.
- 345 Stapley E.O., US patent 4,540,575, 1984.

References *continued*

- 346 Celmer W.D. et al, US patent 4,436,747, 1984.
- 347 Martui M. et al, Agric. Biol. Chem., 1985, 49, 3399.
- 348 Misono H. et al, Agric. Biol. Chem., 1985, 49, 117.
- 349 Ezure Y. et al, Agric. Biol. Chem., 1985, 49 (4) 1119-1125.
- 350 Rose K.E. et al, Appl. and Environ. Microbiol., 1985, 49, 1335.
- 351 Ogata S. et al, Appl. and Environ. Microbiol., 1985, 49, 201.
- 352 Brown R.G. and Embil J.M., Can. J. Microbiol., 1985, 31, 513.
- 353 Brana A.F. et al, Can. J. Microbiol., 1985, 31, 736.
- 354 Carly C.E. et al, Dev. in Ind. Microbiol., 1985, 26, 763.
- 355 Monaghan R.L. and Currie S.A., Dev. in Ind. Microbiol., 1985, 26, 787.
- 356 Reuveny S. et al, Develop. Biol. Std., 1985, 60.
- 357 Horan A.C. et al, EP 0.111.531, 1985.
- 358 EP 0.146.523, 1985.
- 359 Ruddock J.C. et al, EP 0.304.792, 1985.
- 360 Hochstein L.L. and Tomlinson G.A., FEMS Microbiol. Lett., 1985, 27, 329.
- 361 Cryz S.J. Jr. et al, Infect. and Immun., 1985, 50, 225.
- 362 Wright A.C. et al, Infect. and Imm., 1985, 50, 922.
- 363 Goetz M.A. et al, J. Antibiot., 1985, 38 (2) 161-168.
- 364 Boeck L.D. et al, J. Antibiot., 1985, 38, 1.
- 365 Nisbet L.J. et al, J. Antibiot., 1985, 38, 133.
- 366 Tunac J.B. et al, J. Antibiot., 1985, 38, 1337.
- 367 Tunac J.M. and Underhill M., J. Antibiot., 1985, 38, 1344.
- 368 Carler G.T. et al, J. Antibiot., 1985, 38, 242.
- 369 Tunac J.B. et al, J. Antibiot., 1985, 38, 460.
- 370 Kubo K. et al, J. Antibiot., 1985, 38, 904.
- 371 Nakatsukasa W.M. et al, J. Antibiot., 1985, 38, 997.
- 372 Sullivan M.A. and Bock R.M., J. Bacteriol., 1985, 161, 377.
- 373 Phadnis S.H. and Berg D.E., J. Bacteriol., 1985, 161, 928.
- 374 Callahan A. M. and Parkinson J.S., J. Bacteriol., 1985, 161, 96.
- 375 Hansen E.B. et al, J. Bacteriol., 1985, 162, 1227.
- 376 Pang P.P., J. Bacteriol., 1985, 163, 1007.
- 377 Slocum M.K. and Parkinson J.S., J. Bacteriol., 1985, 163, 586.
- 378 Nishibuchi M. and Kaper J.B., J. Bacteriol., 1985, 162, 558.
- 379 Jones G.H., J. Bacteriol., 1985, 163, 1215.
- 380 Anonymous, The United States Pharmacopeia, 1985, XXI, 912.
- 381 Tunac J.B. et al, US patent 4,495,286, 1985.
- 382 Kucera C.J., US patent 4,506,017, 1985.
- 383 Goodnow R.A. et al, US patent 4,530,832, 1985.
- 384 Cerini C.P., US patent 4,537,769, 1985.
- 385 Monaghan R.L. and Kong Y., US patent 4,540,575, 1985.
- 386 Celmer W.D. et al, US patent 4,543,334, 1985.
- 387 Celmer W.D. et al, US patent 4,552,843, 1985.
- 388 Masurekar P.S., ACS Symp. Ser. 362, 1986, 20, 242.
- 389 Awata M. et al, Agric. Biol. Chem., 1986, 50, 239.
- 390 Mieschendahl M. et al, Bio Technology., 1986, (4), 802.
- 391 Loison G. et al, Bio Technology., 1986, 4, 433.
- 392 Friis N. and Myers-Keith P., Biotech. and Bioeng., 1986, 28, 21.
- 393 Greenstein J.P. and Winitz M., Chemistry of the Amino Acids. 2, Robt. E. Krieger Publ Co., Malabac, Fla., 1804, 1986.
- 394 Labeda D.P., Dev. in Ind. Microbiol., 1986, 28, 115.
- 395 EP 0.178.106, 1986.
- 396 McMullen J.R., EP 0.208.489, 1986.
- 397 Myers-Keith P. and Cain W.J., EP 0.209.355, 1986.
- 398 Anonymous, Fisher Scientific Co. Catalog, 585 Alpha Dr., Pittsburgh, PA. 15238, 1986, 463.
- 399 Tomlinson G.A. et al, Int. J. System. Bacteriol., 1986, 36, 66.
- 400 Tunac J.B. et al, J. Antibiot., 1986, 39, 192.
- 401 Reynolds D. et al, J. Appl. Bacteriol., 1986, 60, 517.
- 402 Lunsford R.D. and Macrina F.L., J. Bacteriol., 1986, 166, 426.
- 403 Velez D. et al., J. Imm. Meth., 1986, 86, 45.
- 404 Reuveny S. et al, J. Imm. Meth., 1986, 86, 61.
- 405 Reuveny S. et al, J. Imm. Meth., 1986, 86, 53.
- 406 Marshall V.P. et al, J. Ind. Microbiol., 1986, 1, 17.
- 407 Johnson K.A., Methods in Enzymol., 1986, 134, 306.
- 408 Despreaux C.W. et al, ppl. and Environ. Microbiol., 1986, 51, 946.
- 409 Lee M.D et al, US patent 4,568,646, 1986.
- 410 Meyers E. et al, US patent 4,601,904, 1986.
- 411 Michel K.H. and Kastner R.E., US patent 4,604,239., 1986.
- 412 Tunac J.B., US patent 4,615,975, 1986.
- 413 Shioi I. et al, US patent 4,618,580, 1986.
- 414 Yamanobe T. et al, Agric. Biol. Chem., 1987, 51, 65.
- 415 Kane J.F. and Bogosian G., Bio PhAnn., 1987, 26.
- 416 Bullock W.O. et al, Bio Techniques., 1987, 5, 376.
- 417 Martin N. et al, Bio Technology., 1987, 5, 838.
- 418 Allen B.A. and Luli G.W., BioPhAnn., 1987, 38
- 419 Brennan A.J. et al, Biotechnol.Tech., 1987, 1, 169.
- 420 EP 0.257.276, 1987.
- 421 Hamill R.I. et al, EP 0.308.219, 1987.
- 422 Grimble G.K. et al, Gastroenterology, 1987, 92, 136.
- 423 Moore L.V.H. et al, Int. J. Syst. Bacteriol., 1987, 36, 271.
- 424 Carter G.T. et al, J. Antibiot., 1987, 40 (2) 233-236
- 425 Wiley P.F. et al, J. Antibiot., 1987, 40, 195.
- 426 Tsai M. et al, J. Bacteriol., 1987, 169, 5556.
- 427 Jones G.H., J. Bacteriol., 1987, 169, 5575.
- 428 Mishra S.K. et al, J. Ind. Microbiol., 1987, 2, 267.
- 429 Bailey F.J et al, J. Ind. Microbiol., 1987, 2, 47.
- 430 Sussman M., Meth. Cell Biol., 1987, 28, 9.
- 431 Saxena P.K. and King J., Plant-Cell-Tissue-Organ-Cultures, 1987, 9 (1) 61-71.
- 432 Witterland W.F. Rijksinstituut voor Volksgezondheid en Milieuhygiene, 1987.
- 433 Wilson K.E. and Zimmerman S.B., US patent 4,656,036, 1987.
- 434 Boeck L.D., US patent 4,683,204, 1987.
- 435 Nakatsukasa W., US patent 4,695,545, 1987.
- 436 Hamill R.L. and Yao R.C., US patent 4.683,201, 1987.
- 437 Okamura T. et al, Agric. Biol. Chem., 1988, 52, 1433.
- 438 Hayashida S. et al, Agric. Biol. Chem., 1988, 52, 2397.
- 439 Reynolds D. et al, Appl. and Environ. Microbiol., 1988, 54 (7) 1761-1765.
- 440 Reynolds D. et al, Appl. and Environ. Microbiol., 1988, 54, 1761.
- 441 Gilliland S.E. and R. C., Appl. and Environ. Microbiol., 1988, 54, 898.
- 442 A. Otero A. et al, J. Appl. Bacteriol., 1988, 64, 117.
- 443 Duchene M. et al., J. Bacteriol., 1988, 170, 155.
- 444 Inbar L. and Lapidat A., J. Bacteriol., 1988, 170, 4055.
- 445 Krause M. and Marahiel M.A., J. Bacteriol., 1988, 170, 4669.
- 446 Deich R.A. et al, J. Bacteriol., 1988, 170, 489.
- 447 Dicker H.J. and French E.J., poster presented at 45th Ann. Mtg. SIM., 1988,
- 448 Tunac J.B., US patent 4,734,370, 1988.
- 449 Keith P.M. and Cain W.J., US patent 4,762,784, 1988.
- 450 Lee T.M. et al, US patent 4,774,184, 1988.
- 451 Oka M. and Koniski M., US patent 4,777,160, 1988.
- 452 McMullen J.R., US patent 4,788,144, 1988.
- 453 Wang G. et al, Am. Biotech. Lab., 1989, 7, 14.
- 454 Gomez-Lucia E. et al, Appl. and Environ. Microbiol., 1989, 55, 1447.
- 455 Liszewski M.K. et al, Bio Techniques., 1989, 7, 1079.
- 456 Chisholm D., Bio Techniques., 1989, 7, 21.
- 457 Marshall V.P. et al, J. Antibiot., 1989, 42, 132.
- 458 Dolinger D.L. et al, J. Bacteriol., 1989, 171, 4355.
- 459 Martin M. et al, J. Bacteriol., 1989, 171, 2406.
- 460 Beach M.J. and Rodwell V.W., J. Bacteriol., 1989, 171, 2994.
- 461 Duchene M. et al, J. Bacteriol., 1989, 171, 4130.
- 462 Liesack W. and Stackebrandt E., J. Bacteriol., 1989, 171, 5025.
- 463 Cameron B. et al, J. Bacteriol., 1989, 171, 547.
- 464 McCarter L. and Silverman M., J. Bacteriol., 1989, 171, 731.
- 465 Horinouchi S. et al, J. Bacteriol., 1989, 171, 1206.
- 466 Chen T.S et al, J. Ind. Microbiol., 1989, 4, 231.
- 467 Vecht-Lifshitz S.E. and Braun S., Lett. Appl. Microbiol., 1989, 9, 79.

- 468 Kegel M.A. and Wallace D.L., US patent 4,806,479, 1989.
- 469 Eaton T.E. et al, US patent 4,800,157, 1989.
- 470 Brehm R.D. et al, Appl. and Environ. Microbiol., 1990, 56, 1067.
- 471 Durham D.R., Appl. and Environ. Microbiol., 1990, 56, 2277.
- 472 Ford L.M. et al, Appl. and Environ. Microbiol., 1990, 56, 3511.
- 473 Klein R.C. et al, Bio Techniques., 1990, 9, 196.
- 474 Katz E., Dev. in Ind. Microbiol., 1990, 31, 1.
- 475 Lee Wong A.C. and Bergdoll M.S., Infect. and Immun., 1990, 58, 1026.
- 476 Maiiese W.M. et al, J. Antibiot., 1990, 43, 1059.
- 477 "Schaeffer J.M. et al, J. Antibiot., 1990, 43, 1179.
- 478 Sawada Y. et al, J. Antibiot., 1990, 43, 1223.
- 479 Jackson M. et al., J. Antibiot., 1990, 43, 223.
- 480 Naruse N. et al, J. Antibiot., 1990, 43, 267.
- 481 Carler G.T. et al, J. Antibiot., 1990, 43, 504.
- 482 Dirlam J.P. et al, J. Antibiot., 1990, 43, 668.
- 483 Tomita K. et al, J. Antibiot., 1990, 43, 755.
- 484 Sawada Y. et al, J. Antibiot., 1990, 43, 771.
- 485 "Tatsuta, K. et al, J. Antibiot., 1990, 43, 909.
- 486 Oka M. et al, J. Antibiot., 1990, 43, 967.
- 487 Pearce C.J. et al, J. Antibiot., 1990, 44, 1247.
- 488 Lam Y.K.T. et al, J. Antibiot., 1990, 44, 613.
- 489 Wanner B.L. and Boline J.A., J. Bacteriol., 1990, 172, 1186.
- 490 Ogura T. et al, J. Bacteriol., 1990, 172, 1562.
- 491 Dekleva M.L. and Dasgupta B.R., J. Bacteriol., 1990, 172, 2498.
- 492 Showalter R.E et al, J. Bacteriol., 1990, 172, 2946.
- 493 "Fulks, K.A. et al, J. Bacteriol., 1990, 172, 310.
- 494 Matsushita O. et al, J. Bacteriol., 1990, 172, 3620.
- 495 Kariyama R. et al, J. Bacteriol., 1990, 172, 3718.
- 496 Tsaneva I.R. and Weiss B., J. Bacteriol., 1990, 172, 4197.
- 497 Emori M. et al, J. Bacteriol., 1990, 172, 4901.
- 498 Kobayashi M. et al, J. Bacteriol., 1990, 172, 4807.
- 499 Tomcsanyi T. et al, J. Bacteriol., 1990, 172, 6348.
- 500 Kenny G.E. et al., J. Clin. Microbiol., 1990, 28, 2087.
- 501 Aldwin L. et al, J. Ind. Microbiol., 1990, 5, 239.
- 502 Culler W.P. et al, J. Ind. Microbiol., 1990, 5, 365.
- 503 Dirlam J.P. et al, J. Ind. Microbiol., 1990, 6, 135.
- 504 Leung J.C. et al, Presentation at Am. Chem. Soc. Nat. Mtg., Boston, 1990.
- 505 Vandenberghe P.A. and Kunka B.S., US patent 4,942,032, 1990.
- 506 Mizrahi A. and Lazar A., Anim.Cell-Cult.Prod.Biol., 1991, 159.
- 507 EP 0,456,510, 1991.
- 508 Kobayashi M. et al, FEMS-Microbiol.Lett., 1991, 77 (1), 121
- 509 Schouls L.M. et al, Infect. and Imm., 1991, 59, 3536.
- 510 Hufnagle W.O et al, Infect. and Immun., 1991, 59, 2126.
- 511 Regassa L.B. et al, Infect. and Immun., 1991, 59, 955.
- 512 Vaughn J.L., In-Vitro., 1991, 27, 3, Pt.2, 37A.
- 513 Hensens O.D. et al, J. Antibiot., 1991, 44 (2), 249
- 514 Joshua H. et al, J. Antibiot., 1991, 44, 366.
- 515 Lam Y.K.T. et al, J. Antibiot., 1991, 44, 6, 613.
- 516 Smitka T.A. et al, J. Antibiot., 1991, 45, 433.
- 517 Dombrowski A.W. et al, J. Antibiot., 1991, 45, 671.
- 518 Yamauchi N. and Kakinuma K., J. Antibiot., 1991, 45, 767.
- 519 Kelly K.S. et al, J. Bacteriol., 1991, 173, 2297.
- 520 Kulakauskas S. et al, J. Bacteriol., 1991, 173, 2633.
- 521 Schauer A.T. et al, J. Bacteriol., 1991, 173, 5060.
- 522 Schulman C.A. et al, J. Biotech., 1991, 21, 109.
- 523 Cheng P. et al, J. Ind. Microbiol., 1991, 7, 27.
- 524 de Roubin M.R. et al, J. Ind. Microbiol., 1991, 8, 147.
- 525 Marcin C. et al, J. Ind. Microbiol., 1991, 8, 157.
- 526 Hegde V.R. et al, J. Ind. Microbiol., 1991, 8, 187.
- 527 Hensens O.D. et al, J. Org.Chem., 1991, 56, 10, 3399.
- 528 Mazars C. et al, Prod.Biol.Anim.Cells-Cult., 1991, ESACT 10 Meet., 85.
- 529 US Patent 4,997,849, 1991.
- 530 Christensen B.E. et al, WO 9105858, 1991.
- 531 Baker P.J., WO 9106552, 1991.
- 532 Schmid G. et al, Anim.Cell-Technol., 1992, 203.
- 533 Schmid G. and Keller T., Anim.Cell-Technol., 1992, 285.
- 534 Whelan S.M. et al, Appl. and Environ. Microbiol., 1992, 58, 2345.
- 535 Marugg J.D. et al, Appl. and Environ. Microbiol., 1992, 58, 2360.
- 536 Martens D.E. et al, Biotechnol.Bioeng., 1992, 39, 9, 891.
- 537 EP 0,467,715, 1992.
- 538 EP 0,472,011, 1992.
- 539 EP 0,485,091, 1992.
- 540 Chen T.S. et al, J. Antibiot., 1992, 45 (1) 118.
- 541 O'Sullivan J. et al, J. Antibiot., 1992, 45 (3) 306-312
- 542 Schwartz-R.E. et al, J. Antibiot., 1992, 45, 12, 1853-66.
- 543 Boeck L.D. et al, J. Antibiot., 1992, 45, 8, 1222-30.
- 544 Steenbergen S.M. et al, J. Bacteriol., 1992, 174, 1099.
- 545 Regassa L.B. and Betley M.J., J. Bacteriol., 1992, 174, 5095.
- 546 Eident-Wilkinson B. et al, J. Bacteriol., 1992, 174, 477.
- 547 Regassa L.B. et al, Inf. and Immun., 1992, 60, 3381.
- 548 Schots A. et al, Techniques in Fish Immununology, Fair Haven, N.J., 1992.
- 549 US patent 5,132,216, 1992.
- 550 Uyeda K.A. et al, US patent 6,294,184, 1992.
- 551 Monget D. et al, WO 9212259, 1992.
- 552 Sutar I.I. et al, World-J.Microbiol.Biotechnol., 1992, 8, 3, 254.
- 553 Babb B.L. et al, FEMS Microbiol. Lett., 1993, 114, 3, 343.
- 554 Tkacz J.S. et al, Ind.Microbiol., 1993, 11, 2, 95-103.
- 555 Hayashi K. et al, J. Antibiot., 1993, 46, 12, 1904.
- 556 Chen T.S. et al, J. Nat.Prod., 1993, 56, 5, 755.
- 557 Parrado J. et al, Process-Biochem., 1993, 28, 2, 109.
- 558 US patent 5,185,251, 1993.
- 559 US patent 5,215,900, 1993.
- 560 US patent 5,217,882, 1993.
- 561 US patent 5,218,215, 1993.
- 562 US patent 5,221,625, 1993.
- 563 GB 2,280,435, 1994.
- 564 Hayashi K. et al, J. Antibiot., 1994, 47 (10) 1104.
- 565 Barrow C.J. et al, J. Antibiot., 1994, 47 (11) 1182-1187
- 566 Carter G.T. et al, J. Antibiot., 1994, 47 (12) 1549-1553
- 567 Chen T.S. et al, J. Antibiot., 1994, 47 (12) 1557-1559
- 568 Bernan V.S. et al, J. Antibiot., 1994, 47 (12) 434-1441
- 569 Lam Y.K.T. et al, J. Antibiot., 1994, 47, 6, 724.
- 570 Steinberg D.A. et al, J. Antibiot., 1994, 47, 8, 887.
- 571 JP 06,345,796, 1994.
- 572 Baumert A. et al, Planta-Med., 1994, 60, 2, 143.
- 573 US patent 5,376,367, 1994.
- 574 WO 9418207, 1994.
- 575 Bonarius H.P.J. et al, Biotechnol.Bioeng., 1995, 45, 6, 524.
- 576 LaPorte T.L. et al, Biotechnol.Tech., 1995, 9, 11, 837.
- 577 Gerson D.F. and Xiao X., CA 2,129,416, 1995.
- 578 Russotti G. et al, J. Chem.Technol.Biotechnol., 1995, 63, 1, 37.
- 579 US patent 5,426,108, 1995.
- 580 WO 9529231, 1995.
- 581 Donaldson M.S. and Shuler M.L., Abstr.Pap.Am.Chem.Soc., 1996, 211 Meet., Pt.1, BIOT057.
- 582 Chen R.H. et al, Antibiot., 1996, 49, 6, 596.
- 583 Bonarius H.P.J. et al, Biotechnol.Bioeng., 1996, 50, 3, 299-318.
- 584 Son K.H. et al, Biotechnol.Lett., 1996, 18, 11, 1297.
- 585 Pasupuleti V.K. and Mohandas D., In-Vitro., 1996, 32, 9, 584
- 586 Roggo B.E. et al, J. Antibiot., 1996, 49 (1) 13-19.
- 587 Roberge C. et al, J. Ferment.Bioeng., 1996, 81, 6, 530.
- 588 Schlaeger E.J., J. Immunol.Methods., 1996, 194, 2, 191-199.
- 589 Kamigiri K. et al, J.Antibiot., 1996, 49, 2, 136.
- 590 US patent 5,536,660, 1996.
- 591 WO 9603428, 1996.

References *continued*

- 592 Lambert N. and Merten O.W., *Biotechnol.Bioeng.*, 1997, 2, 54, 165-80.
- 593 Gu X. et al, *Biotechnol.Bioeng.*, 1997, 56, 4, 353-360.
- 594 Shaffiee A. et al, *J. Antibiot.*, 1997, 50, 5, 418-23.
- 595 Yoo I.K. et al, *J. Ferment.Bioeng.*, 1997, 84 (2) 172-175.
- 596 Ishida Y. et al, *J. Ferment.Bioeng.*, 1997, 84, 4, 348-350.
- 597 JP 09.028.373, 1997.
- 598 JP 09.052.842, 1997.
- 599 US patent 5.496.931, 1997.
- 600 Su Y.C. and Wong A.C.L., US patent 5.601.830, 1997.
- 601 Keen M.J. and Rapson N.T., US patent 5.633.162, 1997.
- 602 Primrose S. et al, US patent 5.691.173, 1997.
- 603 Wan N.C. and Goodrick J.C., US patent 5.691.202, 1997.
- 604 Kapat A. et al, *Biotechnol.Lett.*, 1998, 20, 3, 319-23.
- 605 Donaldson M.S. and Shuler M.L., *Biotechnol.Prog.*, 1998, 14 (4), 573-579.
- 606 Dekker K.A. et al, *J. Antibiot.*, 1998, 51, 2, 145-52.
- 607 JP 10.042.900, 1998.
- 608 Etcheverry T. and Ryll T., US patent 5.705.364, 1998.
- 609 US patent 5.710.251, 1998.
- 610 Blom W. et al, US patent 5.741.705, 1998.
- 611 US patent 5.750.695, 1998.
- 612 Kop S.Y. et al, US patent 5.767.069 / 6.255.100, 1998.
- 613 Denoya C.D. et al, US patent 5.776.735, 1998.
- 614 Cote G.L. and Wyckoff H., US patent 5.889.179, 1998.
- 615 Byrne K.M. et al, US patent 5.945.317, 1998.
- 616 Bulthuis B.A. et al, US patent 6.358.716, 1998.
- 617 Gorfien S.F. et al, WO 9808934, 1998.
- 618 Olivier R. et al, WO 9854296, 1998.
- 619 Gupta P., WO 9854352, 1998.
- 620 Platas G. et al, *World-J.Microbiol.Biotechnol.*, 1998, 14 (4), 521-527.
- 621 Nyberg G.B. et al, *Biotechnol.Bioeng.*, 1999, 62 (3), 324-35"
- 622 Brockhaus M. et al, EP 0.939.121, 1999.
- 623 Briggs B.S. et al, *J. Ind.Microbiol.Biotechnol.*, 1999, 23, 3.
- 624 Damhus et al, US patent 5.855.621, 1999.
- 625 US patent 5.856.142, 1999.
- 626 Leqoux R. et al, US patent 5.856.142, 1999.
- 627 Jansen K.-U. et al, US patent 5.888.516, 1999.
- 628 Chartrain M.M., US patent 5.900.368, 1999.
- 629 US patent 5.919.671, 1999.
- 630 Rosson R.A. et al, WO 9932604, 1999.
- 631 Franek F. et al, *Biotechnology Progress*, 2000, 16 (5), 688-692.
- 632 Heideman R. et al, *Cytotechnology*, 2000, 32 (2), 157-167.
- 633 Sugie Y. et al, US 6,277,860, 2000.
- 634 Perkins J.B. et al, EP 1.001.026, 2000.
- 635 Chou C. and Hou J.W., *Int. J. Food Microbiol.*, 2000, 56 (2,3) 113-121.
- 636 Kanazaki H. et al, *J. Antibiot.*, 2000, 53 (1), 58-62.
- 637 Kwon J.Y. et al, *J. Antibiot.*, 2000, 53 (8), 799-806.
- 638 Igarashi Y. et al, *J. Antibiot.*, 2000, 53 (9), 928-933.
- 639 Kirk J.M. et al, *J. Environmental Science and Health*, 2000, 35 (1), 103-120.
- 640 Ohtani T. et al, JP 2000,072765, 2000.
- 641 Koh W.L. and Loh C.S., *Plant Cell Rep.*, 2000, 19 (12) 1177-1183.
- 642 US patent 6.087.126, 2000.
- 643 Price P.J. et al, US patent 6.103.529, 2000.
- 644 Hofmann K.J. et al, US patent 6.159.729, 2000.
- 645 Leet J.E. et al, WO 0003722, 2000.
- 646 Junker B. et al, *Bioprocess and Biosystem Engineering*, 2001, 24 (1) 39-49.
- 647 Chen J., CN 1.270.813, 2001.
- 648 Speck D., EP 1.074.612, 2001,
- 649 Ikonomou L. et al, *In Vitro Cell. Dev. Biol. Anim.*, 2001, 37 (9) 549-559
- 650 Zhang J. and Greasham R.L., US patent 6.232.111, 2001.
- 651 Leung W.S. et al, US patent 6.258.560, 2001.
- 652 Vandenberghe P.A. et al, US patent 6.277.374, 2001.
- 653 Mukeji P.G. et al, US patent 6.287.866, 2001.
- 654 Tuite M.F. et al, US patent 6.291.205, 2001.
- 655 Hohmann H.P. et al, US 6,395,710, 2001.
- 656 Charles S. et al, WO 0134191, 2001.
- 657 Reiter et al, WO 01/23527, 2001.
- 658 Demain A.L. and Fang A., WO200105997 / US patent 6.558.926, 2001.
- 659 Ikonomou L. et al, *Biotechnol.Lett.*, 2002, 24 (12) 965-969
- 660 Schlaeger E.J. et al, EP 1.182.208, 2002.
- 661 Keen M.J. et al, EP 1.221.476, 2002,
- 662 Sakemi S. et al, *J. Antibiot.*, 2002, 55 (1) 6-18
- 663 Lam K. et al, *J. Ind.Microbiol.Biotechnol.*, 2002, 29 (6) 299-302
- 664 Morris A.E. and Reddy P., US 6,413,744, 2002.
- 665 Chu D. et al, US patent 2002/011993,8, 2002.
- 666 Mehta N.M. et al, US 6,210,925, 2002.
- 667 Terman D.S., US 5,728,388, 2002.
- 668 Hoeffkes H. et al, US 9,864,616, 2002.
- 669 Chiarello R.H. et al, US 9,887,149, 2002.
- 670 Hirayama R.C. et al, US patent 6.350.604, 2002.
- 671 Shibuya K. et al, US patent 6.406.909, 2002.
- 672 Tsuchiya R., US patent 6.413.501, 2002.
- 673 Katinger H. et al, WO 02/24876, 2002.
- 674 Wilson G. et al, WO 02/29084, 2002.
- 675 Wan N. et al, WO 02/40686, 2002.
- 676 Baca M. et al, US 6,884,879, 2003.
- 677 Marteijn R.C. et al, *Biotech. and Bioeng.*, 2003, 81(3) 269-278
- 678 Tarnowski J. et al, WO2004085474, 2004.
- 679 Kwon M., Dojima T. and Park E. et al, *Biotechnol.and Applied Biochem.*, 2005, 42 (1) 1-7.
- 680 Baer et al, US patent 4,014,991, 1977.
- 681 US patent 4877615, 1989.
- 682 Fratelli F. et al, *Biotechnol. Prog.*, 2005, 21, 756-761.
- 683 G. Parc et al, *Phytochemistry* 59, 2002, 725-730
- 684 Reeves P. et al, *Proc. Natl. Acad. Sci USA*, 2002, 99(21) 13413-13418.
- 685 O.L. Gamborg, R.A. Miller and K.Ojima, *Exper Cell Research* 50, 1968, 151-158.
- 686 G. Parc et al, *Plant Cell Rep* 26, 145-152.
- 687 Fang et al, *Letters in Applied Microbiol.*, 2006, (43),360-363.
- 688 Demain et al., *Letters in Applied Microbiol.*, 2007, (45),635-63
- 689 WO 042542, 2006.
- 690 Demain et al., *Vaccine* 23 2005, 5420-5423.

Literature references cited

- Barr, K.A., Hopkins, S.A. and Sreekrishna, K. (1992) Protocol for efficient secretion of HAS developed from *Pichia pastoris*. Pharm. Eng. 12: 45-51.
- Blattner, F.R. et al. (1977) Science 196, 161.
- Burteau CC, Verhoeve FR, Mols JE, Ballez J-S, Agathos SN, and Schneider Y-J. (2004) Fortification of a protein-free cell culture medium with plant peptones improves cultivation and productivity of an interferon- γ -producing CHO cell line. In Vitro Cell Dev Biol. Animal. 39: 291-296.
- Fornwald, J.A., Donovan, M.J., Gerber, R, Keller, J., Taylor D.P., Arcuri, E.J. and Brawner, M.E. (1993). Soluble forms of the human T cell receptor CD4 are efficiently expressed by *Streptomyces lividans*. Biotechnology 11: 1031-1036.
- Franek, F and Katinger, H. (2002). Specific effects of synthetic oligopeptides on cultured animal cells. Biotech. Prog. 18, 155-158.
- Galindo E., Bolivar F. & Quintero R. (1990). Maximizing the expression of recombinant proteins in Escherichia coli by manipulation of the culture conditions. J. Ferment. Bioeng. 69:159-165.
- Nagodawithana T.W. & Reed G. (1998) ISBN 0-9646172-8-5. Nutritional requirements of commercially important microorganisms
- Nancib, N., Branolant, H. And Bourdant, C. (1991) Metabolic roles of peptone and yeast extract for the culture of a recombinant strain of *E. coli*. J. Ind. Microbiol. 8: 165-170.
- Siegel R. and Ryu D.D.Y. (1985). Kinetic study of instability of recombinant plasmid pPLc23trpA1 in *E. coli* using two-stage continuous culture system. Biotechnol. Bioeng. 27: 28-33.
- Schlaeger E.J. (1996). The protein hydrolysate Primatone™ RL is a cost effective multiple growth promoter of mammalian cell culture in serum containing and serum free media and displays anti-apoptosis properties. J. Immunol. Methods, 194, 191-199.
- Wang Z. and Da Silva N.A. (1993) Improved protein synthesis and secretion through medium enrichment in a stable recombinant yeast strain. Biotechnol. Bioeng. 42: 95-102.

Amicase, Amisoy, Edamin, Hy-Case, Hy-Express, Hy-Pea, Hy-Soy, Hy-Rice, HyPep, Hy-Wheat, Hy-Yest, Kerry, N-Z-Amine, N-Z-Case, N-Z-Soy, Pepticase, Primatone, Sheffield, Sheftone and UltraPep are trademarks of Kerry Group.

Patents:

EP 0 811 056 B1 - Peptones for tissue culture media

US 5,741,705 - Method for in vitro growth of eucaryotic cells using low molecular weight peptones

US 5,885,835 - Kit for in vitro cell growth of eucaryotes using low molecular wieght peptone containing L-glutamine

Please note that the information herein is believed to be accurate and complete. However, no warranty or guarantee whatsoever is made or is to be implied with respect to such information or with respect to any other product, method, performance or apparatus referred herein. Nothing herein shall be construed as a recommendation to use, manufacture, sell or act in violation of any intellectual property rights or legal requirements.

KERRY



kerry.com/pharma

© Kerry 2021